



DICOM Conformance Statement

Vitrexa Connection 9.2

July 10, 2024

Canon Medical Informatics Incorporated
19 Regina Street North
Waterloo, Ontario
N2J 2Z9
Canada

Canon Medical shall not be liable for errors contained herein or for incidental or consequential damages in connection with the furnishing, performance or use of this publication. Canon Medical reserves the right to revise this publication and to make changes to its content at any time, without obligation to notify any person or entity of such revisions and changes. This publication may only be used in connection with the promotion, sales, installation and use of Canon Medical products.

mi.medical.canon

Preface

Vitre Connection is a self-contained networked information system used for receiving, archiving and displaying of diagnostic medical images. The system conforms to the DICOM 3.0 standard to allow the sharing of medical information with other digital imaging systems.

NETWORK SERVICES

DICOM SOP Class Name	User of Service (SCU)	Provider of Service (SCP)	Display
Transfer			
12-lead ECG Waveform Storage	Yes	Yes	No
Acquisition Context SR Storage	Yes	Yes	No
Advanced Blending Presentation State Storage	Yes	Yes	No
Ambulatory ECG Waveform Storage	Yes	Yes	No
Arterial Pulse Waveform Storage	Yes	Yes	No
Autorefraction Measurements Storage	Yes	Yes	No
Basic Structured Display Storage	Yes	Yes	No
Basic Text SR	Yes	Yes	No
Basic Voice Audio Waveform Storage	Yes	Yes	No
Blending Softcopy Presentation State Storage	Yes	Yes	No
Breast Tomosynthesis Image Storage	Yes	Yes	No
Breast Projection X-Ray Image Storage – For Presentation	Yes	Yes	No
Breast Projection X-Ray Image Storage – For Processing	Yes	Yes	No
Cardiac Electrophysiology Waveform Storage	Yes	Yes	No
Chest CAD SR	Yes	Yes	No
Colon CAD SR	Yes	Yes	No
Color Softcopy Presentation State Storage	Yes	Yes	No
Corneal Topography Map Storage	Yes	Yes	No
Compositing Planar MPR Volumetric Presentation State Storage	Yes	Yes	No
Comprehensive SR	Yes	Yes	No
Comprehensive 3D SR	Yes	Yes	No
Computed Radiography Image Storage	Yes	Yes	No

Content Assessment Results Storage	Yes	Yes	No
CT Image Storage	Yes	Yes	No
CT Performed Procedure Protocol Storage	Yes	Yes	No
Deformable Spatial Registration Storage	Yes	Yes	No
Digital Intra-oral X-Ray Image Storage - For Presentation	Yes	Yes	No
Digital Intra-oral X-Ray Image Storage - For Processing	Yes	Yes	No
Digital Mammography X-Ray Image Storage - For Presentation	Yes	Yes	No
Digital Mammography X-Ray Image Storage - For Processing	Yes	Yes	No
Digital X-Ray Image Storage - For Presentation	Yes	Yes	No
Digital X-Ray Image Storage - For Processing	Yes	Yes	No
Encapsulated CDA Storage	Yes	Yes	No
Encapsulated PDF Storage	Yes	Yes	No
Encapsulated STL Storage	Yes	Yes	No
Enhanced CT Image Storage	Yes	Yes	No
Enhanced MR Color Image Storage	Yes	Yes	No
Enhanced MR Image Storage	Yes	Yes	No
Enhanced PET Image Storage	Yes	Yes	No
Enhanced SR	Yes	Yes	No
Enhanced US Volume Storage	Yes	Yes	No
Enhanced XA Image Storage	Yes	Yes	No
Enhanced XRF Image Storage	Yes	Yes	No
Extensible SR Storage	Yes	Yes	No
General Audio Waveform Storage	Yes	Yes	No
General ECG Waveform Storage	Yes	Yes	No
Grayscale Softcopy Presentation State Storage	Yes	Yes	No
Grayscale Planer MPR Volumetric Presentation State Storage	Yes	Yes	No
Hanging Protocol Storage	Yes	Yes	No
<i>Hardcopy Color Image Storage (Retired)</i>	Yes	Yes	No
<i>Hardcopy Grayscale Image Storage (Retired)</i>	Yes	Yes	No
Hemodynamic Waveform Storage	Yes	Yes	No

Implantation Plan SR Document Storage	Yes	Yes	No
Intraocular Lens Calculation Storage	Yes	Yes	No
Intravascular Optical Coherence Tomography Image Storage – For Presentation	Yes	Yes	No
Intravascular Optical Coherence Tomography Image Storage – For Processing	Yes	Yes	No
Key Object Selection	Yes	Yes	No
Keratometry Measurements Storage	Yes	Yes	No
Legacy Converted Enhanced CT Image Storage	Yes	Yes	No
Legacy Converted Enhanced MR Image Storage	Yes	Yes	No
Legacy Converted Enhanced PET Image Storage	Yes	Yes	No
Lensometry Measurements Storage	Yes	Yes	No
Macular Grid Thickness and Volume Report	Yes	Yes	No
Mammography CAD SR	Yes	Yes	No
MR Image Storage	Yes	Yes	No
MR Spectroscopy Storage	Yes	Yes	No
Multi-frame Grayscale Byte Secondary Capture Image Storage	Yes	Yes	No
Multi-frame Grayscale Word Secondary Capture Image Storage	Yes	Yes	No
Multi-frame Single Bit Secondary Capture Image Storage	Yes	Yes	No
Multi-frame True Color Secondary Capture Image Storage	Yes	Yes	No
Multiple Volume Rendering Volumetric Presentation State Storage	Yes	Yes	No
Nuclear Medicine Image Storage	Yes	Yes	No
<i>Nuclear Medicine Image Storage (Retired)</i>	Yes	Yes	No
Ophthalmic Axial Measurements Storage	Yes	Yes	No
Ophthalmic Optimal Coherence Tomography En Face Image Storage	Yes	Yes	No
Ophthalmic Optimal Coherence Tomography B-scan Volume Analysis Storage	Yes	Yes	No
Ophthalmic Photography 16 Bit Image Storage	Yes	Yes	No
Ophthalmic Photography 8 Bit Image Storage	Yes	Yes	No
Ophthalmic Thickness Map Storage	Yes	Yes	No
Ophthalmic Tomography Image Storage	Yes	Yes	No

Ophthalmic Visual Field Static Perimetry Measurements Storage	Yes	Yes	No
Parametric Image Storage	Yes	Yes	No
Parametric Map Storage	Yes	Yes	No
Patient Radiation Dose SR Storage	Yes	Yes	No
Performed Imaging Agent Administration SR Storage	Yes	Yes	No
Planning Imaging Agent Administration SR Storage	Yes	Yes	No
Positron Emission Tomography Image Storage	Yes	Yes	No
Procedure Log	Yes	Yes	No
Pseudo-Color Softcopy Presentation State Storage	Yes	Yes	No
Radiopharmaceutical Radiation Dose SR	Yes	Yes	No
Raw Data Storage	Yes	Yes	No
Real World Value Mapping Storage	Yes	Yes	No
Respiratory Waveform Storage	Yes	Yes	No
RT Beams Delivery Instruction Storage	Yes	Yes	No
RT Beams Treatment Record Storage	Yes	Yes	No
RT Brachy Application Setup Delivery Instruction Storage	Yes	Yes	No
RT Brachy Treatment Record Storage	Yes	Yes	No
RT Dose Storage	Yes	Yes	No
RT Image Storage	Yes	Yes	No
RT Ion Beams Treatment Record Storage	Yes	Yes	No
RT Ion Plan Storage	Yes	Yes	No
RT Physician Intent Storage	Yes	Yes	No
RT Plan Storage	Yes	Yes	No
RT Segment Annotation Storage	Yes	Yes	No
RT Structure Set Storage	Yes	Yes	No
RT Treatment Summary Record Storage	Yes	Yes	No
Secondary Capture Image Storage	Yes	Yes	No
Segmentation Storage	Yes	Yes	No
Segmented Volume Rendering Volumetric Presentation State Storage	Yes	Yes	No
Simplified Adult Echo SR Storage	Yes	Yes	No

Spatial Fiducials Storage	Yes	Yes	No
Spatial Registration Storage	Yes	Yes	No
Spectacle Prescription Report Storage	Yes	Yes	No
Stereometric Relationship Storage	Yes	Yes	No
Standalone Overlay Storage (Retired)	Yes	Yes	No
Standalone Curve Storage (Retired)	Yes	Yes	No
Standalone Modality LUT Storage (Retired)	Yes	Yes	No
Standalone VOI LUT Storage (Retired)	Yes	Yes	No
Standalone PET Curve Storage (Retired)	Yes	Yes	No
Subjective Refraction Measurements Storage	Yes	Yes	No
Surface Scan Mesh Storage	Yes	Yes	No
Surface Scan Point Cloud Storage	Yes	Yes	No
Surface Segmentation Storage	Yes	Yes	No
Tractography Results Storage	Yes	Yes	No
Ultrasound Image Storage	Yes	Yes	No
<i>Ultrasound Image Storage (Retired)</i>	Yes	Yes	No
Ultrasound Multi-frame Storage	Yes	Yes	No
<i>Ultrasound Multi-frame Storage (Retired)</i>	Yes	Yes	No
Video Endoscopic Image Storage	Yes	Yes	No
Video Microscopic Image Storage	Yes	Yes	No
Video Photographic Image Storage	Yes	Yes	No
Visual Acuity Measurements Storage	Yes	Yes	No
VL Endoscopic Image Storage	Yes	Yes	No
<i>VL Image Storage (Retired)</i>	Yes	Yes	No
VL Microscopic Image Storage	Yes	Yes	No
VL Photographic Image Storage	Yes	Yes	No
VL Slide-Coordinates Microscopic Image Storage	Yes	Yes	No
VL Whole Slide Microscopy Image Storage	Yes	Yes	No
VOI LUT Box (Retired)	Yes	Yes	No
Volume Rendering Volumetric Presentation State Storage	Yes	Yes	No

Wide Field Ophthalmic Photography Stereographic Projection Image Storage	Yes	Yes	No
Wide Field Ophthalmic Photography 3D Coordinates Image Storage	Yes	Yes	No
X-Ray 3D Angiographic Image Storage	Yes	Yes	No
X-Ray 3D Craniofacial Image Storage	Yes	Yes	No
<i>X-Ray Angiographic Bi-Plane Image Storage (Retired)</i>	Yes	Yes	No
X-Ray Angiographic Image Storage	Yes	Yes	No
X-Ray Radiation Dose SR	Yes	Yes	No
X-Ray Radiofluoroscopic Image Storage	Yes	Yes	No
XA/XRF Grayscale Softcopy Presentation State Storage	Yes	Yes	No
Query/Retrieve			
Patient Root Q/R – FIND	Yes	Yes	N/A
Patient Root Q/R - MOVE	Yes	Yes	N/A
Study Root Q/R – FIND	Yes	Yes	N/A
Study Root Q/R - MOVE	Yes	Yes	N/A
Modality Worklist			
Modality Worklist – FIND	Yes	Yes	N/A
Verification			
Verification	Yes	Yes	N/A

Table of Contents

1	Introduction	13
1.1	Revision History	13
1.2	Audience	13
1.3	Remarks	14
1.4	Definitions, Terms and Abbreviations	14
1.5	References	15
2	Networking	16
2.1	Implementation Model	16
2.1.1	Application Data Flow	16
2.1.2	Functional Definitions of AEs	18
2.1.2.1	Functional Definition of Vitrea Connection Storage Server AE	18
2.1.2.2	Functional Definition of Vitrea Connection Query/Retrieve Server AE	18
2.1.2.3	Functional Definition of Modality Worklist AE	19
2.1.2.4	Functional Definition of Proxy Server AE	19
2.1.2.5	STOW-RS Specification	19
2.1.2.6	QIDO-RS Specification	19
2.1.2.7	WADO-URI Specification	19
2.1.2.8	WADO-RS Specification	19
2.2	AE Specifications	19
2.2.1	Vitrea Connection Storage Server AE Specification	26
2.2.1.1	SOP Classes	26
2.2.1.2	Association Establishment Policy	26
2.2.1.2.1	General	26
2.2.1.2.2	Number of Associations	26
2.2.1.2.3	Asynchronous Nature	26
2.2.1.2.4	Implementation Identifying Information	27
2.2.1.3	Association Acceptance Policy	27
2.2.1.3.1	Activity - Receive Images	27
	Description and Sequencing of Activities	27
	Accepted Presentation Contexts	27
	SOP Specific Conformance	28
	Specific Conformance for the Verification SOP Class	28

Specific Conformance for the Storage SOP Classes	28
Storage Server AE C-STORE Response	29
Specific Conformance for the Storage Commitment SOP Class	30
Storage Server AE N-ACTION Response	30
Storage Server AE N-EVENT Response	31
2.2.1.4 Association Initiation Policy	32
2.2.1.4.1 Activity – Synchronize Objects with an External Peer AE	32
Description and Sequencing of Activities	32
Proposed Presentation Contexts	32
SOP Specific Conformance	32
Specific Conformance for Verification SOP Class	32
2.2.1.4.2 Activity – Archive Objects to an External Peer AE	32
Description and Sequencing of Activities	32
Proposed Presentation Contexts	33
2.2.2 Vitrea Connection Query/Retrieve Server AE Specification	33
2.2.2.1 SOP Classes	33
2.2.2.2 Association Establishment Policy	33
2.2.2.2.1 General	33
2.2.2.2.2 Number of Associations	33
2.2.2.2.3 Asynchronous Nature	34
2.2.2.2.4 Implementation Identifying Information	34
2.2.2.3 Association Initiation Policy	34
2.2.2.3.1 Activity - Send Images Requested by an External Peer AE	34
Description and Sequencing of Activities	34
Proposed Presentation Contexts	34
SOP Specific Conformance	35
Specific Conformance for Verification SOP Class	35
Specific Conformance for Image, SR and KIN SOP Classes	35
2.2.2.3.2 Activity - Retrieve Images Requested by an External Peer AE	36
Description and Sequencing of Activities	36
Proposed Presentation Contexts	36
SOP Specific Conformance	37
Specific Conformance for Verification SOP Class	37
2.2.2.4 Association Acceptance Policy	37

2.2.2.4.1	Activity - Handling Query and Retrieval Requests.....	37
	Description and Sequencing of Activities	37
	Accepted Presentation Contexts.....	37
	SOP Specific Conformance.....	38
	Specific Conformance for Verification SOP Class	38
	Specific Conformance for Query SOP Classes	38
	Specific Conformance for Retrieve SOP Classes	42
2.2.3	Vitrea Connection Modality Worklist AE Specification	44
2.2.3.1	SOP Classes.....	44
2.2.3.2	Association Establishment Policy	44
2.2.3.2.1	General	44
2.2.3.2.2	Number of Associations	44
2.2.3.2.3	Asynchronous Nature.....	44
2.2.3.2.4	Implementation Identifying Information.....	44
2.2.3.3	Association Acceptance Policy	45
2.2.3.3.1	Activity - Handling Query and Retrieval Requests.....	45
	Description and Sequencing of Activities	45
	Accepted Presentation Contexts.....	45
	SOP Specific Conformance.....	45
	Specific Conformance for Verification SOP Class	45
	Specific Conformance for Query SOP Classes	46
2.2.4	Vitrea Connection STOW-RS SCP Specification.....	49
2.2.5	Vitrea Connection QIDO-RS Provider Specification	49
2.2.5.1	QIDO-RS Search for Studies	49
2.2.5.2	QIDO-RS Search for Series	51
2.2.5.3	QIDO-RS Search for Instances	52
2.2.5.4	QIDO-RS Connection Policies.....	53
2.2.5.4.1	General	53
2.2.5.4.2	Number of Connections	53
2.2.5.4.3	Response Status.....	53
2.2.5.4.4	Web Service Endpoint URL.....	53
2.2.6	Vitrea Connection WADO-URI Provider Specification	54
2.2.7	Vitrea Connection WADO-RS Provider Specification	55
2.2.7.1	WADO-RS Retrieve Studies.....	55

2.2.7.2	WADO-RS Retrieve Series	55
2.2.7.3	WADO-RS Retrieve Instance	56
2.2.7.4	WADO-RS Retrieve Frames	57
2.2.7.5	WADO-RS Retrieve Bulk Data	57
2.2.7.6	WADO-RS Retrieve Metadata	58
2.2.7.7	WADO-RS Connection Policies	58
2.2.7.7.1	General	58
2.2.7.7.2	Number of Connections	59
2.2.7.7.3	Web Service Endpoint URL	59
2.2.8	Vitrea Connection Proxy Server AE Specification	59
2.2.8.1	SOP Classes	59
2.2.8.2	Association Establishment Policy	59
2.2.8.2.1	General	59
2.2.8.2.2	Number of Associations	59
2.2.8.2.3	Asynchronous Nature	60
2.2.8.2.4	Implementation Identifying Information	60
2.2.8.3	Association Acceptance Policy	60
2.2.8.3.1	Activity – Receive Images	60
	Description and Sequencing of Activities	60
	Accepted Presentation Contexts	61
	SOP Specific Conformance	61
2.2.8.4	Association Initiation Policy	61
2.2.8.4.1	Activity – Synchronize Objects with an External Peer AE	61
	Description and Sequencing of Activities	61
	Proposed Presentation Contexts	61
	SOP Specific Conformance	62
2.3	Network Interfaces	62
2.3.1	Supported Communication Stacks	62
2.3.1.1	TCP/IP Stack	62
2.3.2	Physical Network Interface	62
2.4	Configuration	62
2.4.1	AE Title/Presentation Address Mapping	62
2.4.1.1	Local AE Titles	62
2.4.1.2	Remote AE Titles	63

2.4.2	Configuration Parameters.....	63
3	Media Interchange.....	64
4	Support of Extended Character Sets.....	64
5	Security.....	64
5.1	Security Profiles	64
5.2	Association Level Security.....	64
5.3	Application Level Security.....	64

1 Introduction

1.1 Revision History

Document Version	Date	Author	Description
0.1	March 13, 2020	Beth Whittle, Joshua Thomas	Initial Draft. Document adapted from Rialto 7.5 DICOM Conformance Statement
1.0	March 23, 2020	Paul Steinbach	Updates for Vitrea Connection 8.0, from review feedback.
1.1	March 26, 2020	Paul Steinbach	Corrections after additional review feedback.
1.2	May 27, 2020	Judith Olivares	Updated for the VC 8.1 release.
1.3	October 16, 2020	Joshua Thomas	Reviewed and updated for the VC 8.2 release.
1.4	March 22, 2021	Joshua Thomas	Reviewed and updated for the VC 8.3 release.
1.5	March 25, 2021	Joshua Thomas	Additional updates made based on reviews.
1.6	June 22, 2021	Joshua Thomas	Updated for the VC 8.4 release.
1.7	February 3, 2022	John McMullen	For VC 8.4.2 release, company name and logo updated to be Canon Medical Informatics Incorporated.
1.8	Nov 29, 2022	Grace Huang	For VC 8.5 release, include updates on DICOMweb.
1.9	July 7, 2023	Grace Huang	For VC 9.0 release, added remarks about VC 9.0 no support/conformance for DICOMweb standards.
2.0	Dec 15, 2023	Grace Huang	For VC 9.1 release. Removed VC9.0 remarks on VC 9.0. DICOMweb conformance included in VC 9.1.
2.1	July 10, 2024	Tim Bartram	Updated with new IHE Radiology Technical Framework version and DICOM version

1.2 Audience

This document is intended for hospital staff, health system integrators, software designers or implementers. It is assumed that the reader is familiar with the terminology and concepts that are used in the DICOM 3.0 standard and the appropriate IHE Profiles.

1.3 Remarks

DICOM, by itself, does not guarantee interoperability. However, the Conformance Statement facilitates a first level validation for interoperability between different applications supporting the same DICOM functionality.

This Conformance Statement is not intended to replace validation with other DICOM equipment to ensure proper exchange of information intended.

The scope of this Conformance Statement is to facilitate communication with Vitrea Connection and other vendors' medical equipment. However, it is not guaranteed to ensure the desired interoperability and successful interconnectivity with existing DICOM systems.

The user should be aware of the following important issues:

- The comparison of different conformance statements is the first step towards assessing interconnectivity between Vitrea Connection and other vendors' equipment.
- Test procedures should be defined to validate the desired level of connectivity.
- The DICOM 3.1 standard will evolve to meet the users' future requirements. Canon Medical reserves the right to make changes to its products or to discontinue its delivery.

1.4 Definitions, Terms and Abbreviations

Definitions, terms and abbreviations used in this document are defined within the different parts of the DICOM standard.

Abbreviations and terms are as follows:

AE	Application Entity
AET	Application Entity Title
DICOM	Digital Imaging and Communications in Medicine
DICOMweb	DICOM Standard for Web-based Medical Imaging
DIMSE	DICOM Message Service Element
GSDF	Grayscale Standard Display Function
GUI	Graphical User Interface
HIS/RIS	Hospital Information System / Radiology Information System
HL7	Health Level 7
IHE	Integrating the Healthcare Enterprise
IHE-TF	Integrating the Healthcare Enterprise Technical Framework
IOD	Information Object Definition
ISO	International Standard Organization
PDU	DICOM Protocol Data Unit
LUT	Look-up Table
P-LUT	Presentation Look-up Table
QIDO	Query based on ID for DICOM Objects
SCP	Service Class Provider
SCU	Service Class User
SOP	DICOM Service-Object Pair
STOW	STore Over the Web
TCP/IP	Transmission Control Protocol/Internet Protocol
TLS	Transport Layer Security
UID	Unique Identifier
VR	Value Representation
WADO	Web Access to DICOM Object
XDS	Cross-Enterprise Document Sharing

1.5 References

- [DICOM] Digital Imaging and Communications in Medicine (DICOM), NEMA PS 3.1 - 3.22, 2024c
- [IHE-TF] IHE Radiology Technical Framework, HIMSS/RSNA, Vol. I – IV, Revision 22.0, June 25, 2024
- [IHE-TF] IHE IT Infrastructure Technical Framework, HIMSS/RSNA, Vol. I – IV, Revision 20, August 4, 2023


2 Networking

2.1 Implementation Model


2.1.1 Application Data Flow

Vitreia Connection is logically divided into different DICOM Application Entities: Storage Server, Query/Retrieve Server, STOW-RS, WADO-URI Provider, and MWL Server. For easy reading, flows for standard DIMSE-C services have been separated from web service entry points.

Application Data Flow Diagram – DIMSE-C

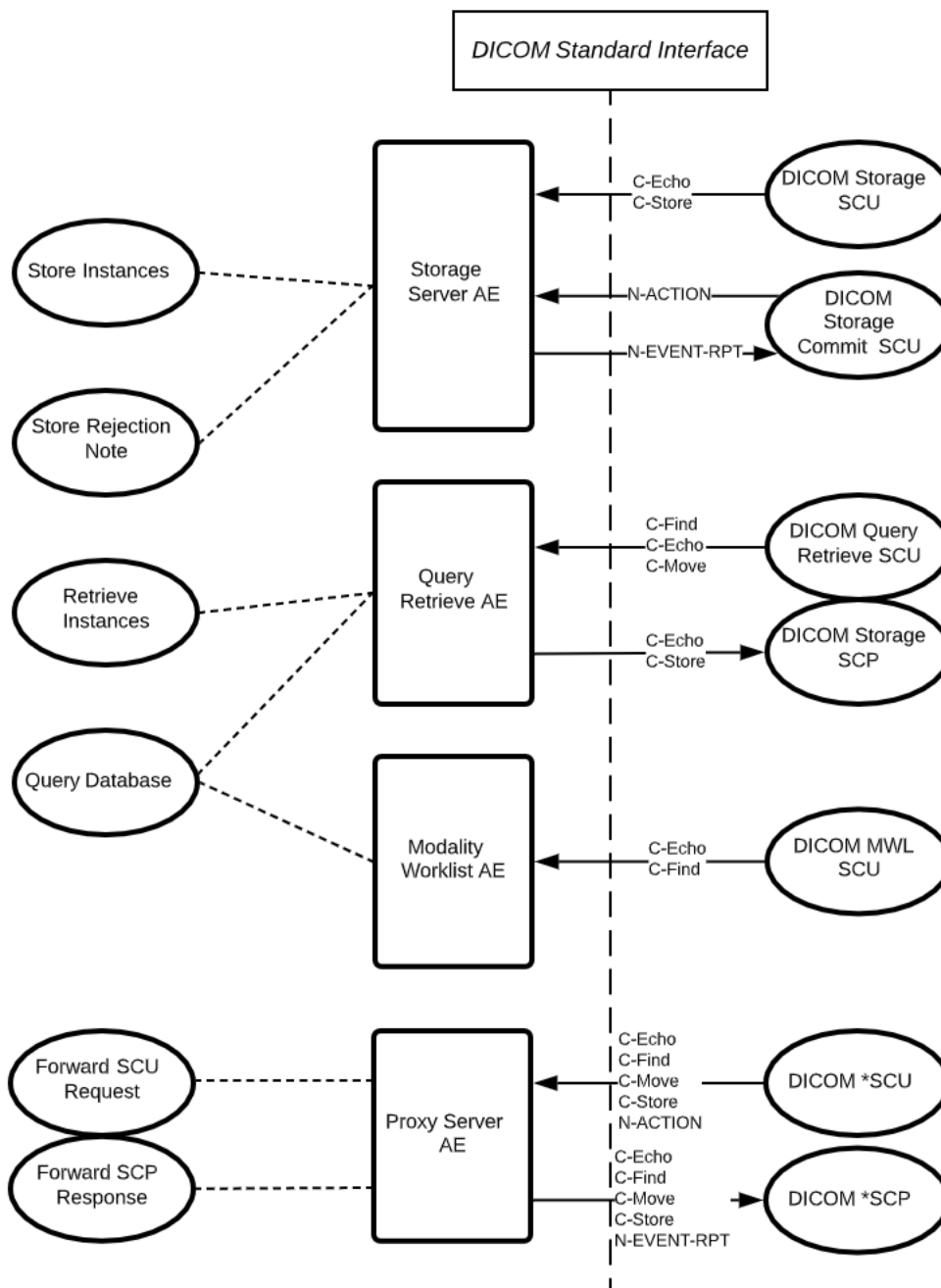

Application Data Flow Diagram – Web Services

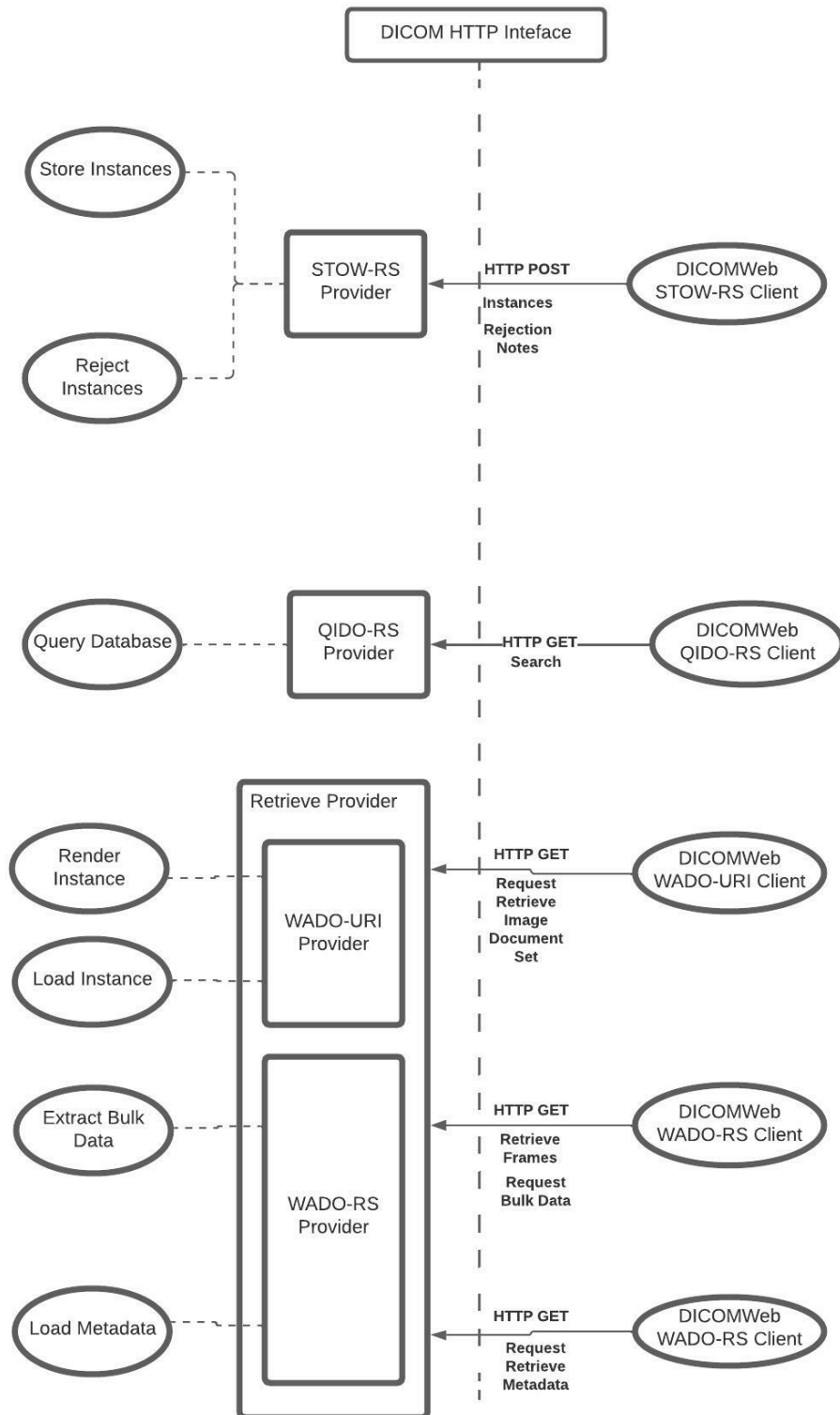

The Storage Server AE implements the DICOM Storage, Storage Commit and Verification Service Classes.

The Query/Retrieve Server AE implements the DICOM Storage, Retrieve and Verification Service Classes.

The Modality Worklist AE implements the DICOM Modality Worklist and Verification Service Classes.

The Proxy Server AE implements the DICOM Storage, Storage Commit, Retrieve and Verification Service Classes.

The STOW-RS SCP AE implements the Store Instances transaction for DICOM Web Services.

The WADO-URI Provider implements the Retrieve Imaging Document Set and Retrieve Rendered Imaging Document Set transactions for DICOM Web Services.

The QIDO-RS Provider implements the Query DB transaction for DICOM Web Services.

The WADO-RS Provider implements Retrieve Frames, Request Bulk Data, and Retrieve Metadata for DICOM Web services.

2.1.2 Functional Definitions of AEs

2.1.2.1 Functional Definition of Vitrea Connection Storage Server AE

Vitrea Connection Storage Server Application Entity waits for another application to connect at the presentation address configured for its Application Entity Title. When another application connects, the Storage Server AE expects it to be a DICOM application.

The Storage Server AE will accept Associations with Presentation Contexts for SOP Classes of the Verification, Storage, and Storage Commitment Service Classes.

Any objects received on such Presentation Contexts will be added to the Object Repository.

As an option, the Storage Server AE can be configured to modify the contents of objects with information found in scheduled procedure steps (i.e., modality worklist reconciliation).

As an option, the Storage Server AE can be configured to extract measurements that are stored as XML strings in objects and to write those measurements to files.

As an option, the Storage Server AE will create a Key Object Selection document and publish it to an XDS Document Repository, thus creating or updating entries in an XDS Registry.

As part of regular maintenance activities, more specifically the synchronization of information that has been received by the Storage Server AE, the Storage Server AE will initiate a request to query for and retrieve objects from the AE that originally stored the object.

As part of regular maintenance activities, more specifically when Vitrea Connection is configured as a cache, the Storage Server AE will initiate a request to store objects to one or more AEs that are serving as a long-term archive.

2.1.2.2 Functional Definition of Vitrea Connection Query/Retrieve Server AE

Vitrea Connection Query/Retrieve Server Application Entity waits for another application to connect at the presentation address configured for its Application Entity Title. When another application connects, the Retrieve Server AE expects it to be a DICOM application.

The Vitrea Connection Query/Retrieve Server AE will accept Associations with Presentation Contexts for SOP Classes of the Verification and Retrieve Service Classes.

The Vitrea Connection Query/Retrieve Server AE responds to queries based on the records stored in its database.

The Vitrea Connection Query/Retrieve Server AE is used to query and retrieve DICOM objects that were stored to the Storage Server AE.

The Vitrea Connection Query/Retrieve Server AE may be participant in a C-Move operation as a C-Store SCU on a separate association.

2.1.2.3 Functional Definition of Modality Worklist AE

Vitreia Connection MWL Server Application Entity waits for another application to connect at the presentation address configured for its Application Entity Title. When another application connects, the Modality Worklist AE expects it to be a DICOM application.

The Modality Worklist AE will accept Associations with Presentation Contexts for SOP Classes of the Verification and Modality Worklist Service Classes.

The Modality Worklist AE responds to queries based on the records stored in its database.

2.1.2.4 Functional Definition of Proxy Server AE

Vitreia Connection Proxy Server Application Entity waits for another application to connect at the presentation address configured for its Application Entity Title. When another application connects, the Proxy Server AE expects it to be a DICOM application.

The Proxy Server AE will accept Associations with Presentation Contexts for SOP Classes of the Verification, Storage Commit, Storage and Query/Retrieve Service Classes.

The Proxy Server AE responds to queries by proxying the Commit/Query/Retrieve/Storage requests to one or more AEs and providing the results to the original requesting AE.

2.1.2.5 STOW-RS Specification

Vitreia Connection STOW-RS AE waits for another application to connect at the URL configured for its Application Entity Title. When another application connects, the STOW-RS AE expects it to be a web application.

The STOW-RS AE will accept web requests for the storage of instances.

The STOW-RS AE will create a Key Object Selection document and publish it to an XDS Document Repository, thus creating or updating entries in an XDS Registry.

2.1.2.6 QIDO-RS Specification

The QIDO-RS Provider provides access to Study, Series, and Instance data of received objects via HTTP GET requests.

2.1.2.7 WADO-URI Specification

Vitreia Connection WADO Server Application Entity waits for another application to connect at the URL configured for its Application Entity Title. When another application connects, the WADO-URI Provider expects it to be a web application.

The WADO-URI Provider will accept web requests for retrieval of instances.

The WADO-URI Provider responds to queries based on the records stored in its database.

The WADO-URI Provider is used to retrieve DICOM objects that were stored to the Storage Server AE. It expects that the requesting application already knows the SOP Instance UIDs of the DICOM objects that it wishes to retrieve.

2.1.2.8 WADO-RS Specification

The WADO-RS Provider provides access to the metadata, the bulk data or the whole DICOM Objects, of a whole Study or Series or an individual object via HTTP GET requests.

2.2 AE Specifications

The following tables define default values that are referenced by the following sections.

The following Transfer Syntaxes are represented by the term [Default TS List]:

Name List	UID List
Implicit VR Little Endian	1.2.840.10008.1.2
Explicit VR Little Endian	1.2.840.10008.1.2.1

The following Transfer Syntaxes are represented by the term [Extended TS List]:

Name List	UID List
Implicit VR Little Endian	1.2.840.10008.1.2
Explicit VR Little Endian	1.2.840.10008.1.2.1
Explicit VR Big Endian	1.2.840.10008.1.2.2
JPEG Baseline (Process 1)	1.2.840.10008.1.2.4.50
JPEG Extended (Process 2 & 4)	1.2.840.10008.1.2.4.51
JPEG Lossless, Non-Hierarchical (Process 14)	1.2.840.10008.1.2.4.57
JPEG Lossless, Non- Hierarchical, First-Order Prediction (Process 14 [Selection Value 1])	1.2.840.10008.1.2.4.70
JPEG-LS Lossless Image Compression	1.2.840.10008.1.2.4.80
JPEG-LS Lossy (Near-Lossless) Image Compression	1.2.840.10008.1.2.4.81
JPEG 2000 Lossless Image Compression	1.2.840.10008.1.2.4.90
JPEG 2000 Lossy Image Compression	1.2.840.10008.1.2.4.91
RLE Lossless	1.2.840.10008.1.2.5

The following SOP Classes are represented by the term [Default SOP Class List]:

Name List	UID List	Transfer Syntax List
<i>VOI LUT Box (Retired)</i>	1.2.840.10008.5.1.1.22	[Default TS List]
<i>Hardcopy Grayscale Image Storage (Retired)</i>	1.2.840.10008.5.1.1.29	[Extended TS List]
<i>Hardcopy Color Image Storage SOP Class (Retired)</i>	1.2.840.10008.5.1.1.30	[Extended TS List]
Computed Radiography Image Storage	1.2.840.10008.5.1.4.1.1.1	[Extended TS List]
Digital X-Ray Image Storage - For Presentation	1.2.840.10008.5.1.4.1.1.1.1	[Extended TS List]
Digital X-Ray Image Storage - For Processing	1.2.840.10008.5.1.4.1.1.1.1.1	[Extended TS List]
Digital Mammography X-Ray Image Storage - For Presentation	1.2.840.10008.5.1.4.1.1.1.2	[Extended TS List]
Digital Mammography X-Ray Image Storage - For Processing	1.2.840.10008.5.1.4.1.1.1.2.1	[Extended TS List]

Digital Intra-Oral X-Ray Image Storage - For Presentation	1.2.840.10008.5.1.4.1.1.1.3	[Extended TS List]
Digital Intra-Oral X-Ray Image Storage - For Processing	1.2.840.10008.5.1.4.1.1.1.3.1	[Extended TS List]
CT Image Storage	1.2.840.10008.5.1.4.1.1.2	[Extended TS List]
Enhanced CT Image Storage	1.2.840.10008.5.1.4.1.1.2.1	[Extended TS List]
Legacy Converted Enhanced CT Image Storage	1.2.840.10008.5.1.4.1.1.2.2	[Extended TS List]
<i>Ultrasound Multi-frame Image Storage (Retired)</i>	1.2.840.10008.5.1.4.1.1.3	[Extended TS List]
Ultrasound Multi-frame Image Storage	1.2.840.10008.5.1.4.1.1.3.1	[Extended TS List]
MR Image Storage	1.2.840.10008.5.1.4.1.1.4	[Extended TS List]
Enhanced MR Image Storage	1.2.840.10008.5.1.4.1.1.4.1	[Extended TS List]
MR Spectroscopy Storage	1.2.840.10008.5.1.4.1.1.4.2	[Extended TS List]
Enhanced MR Color Image Storage	1.2.840.10008.5.1.4.1.1.4.3	[Extended TS List]
Legacy Converted Enhanced MR Image Storage	1.2.840.10008.5.1.4.1.1.4.4	[Extended TS List]
<i>Nuclear Medicine Image Storage (Retired)</i>	1.2.840.10008.5.1.4.1.1.5	[Extended TS List]
<i>Ultrasound Image Storage (Retired)</i>	1.2.840.10008.5.1.4.1.1.6	[Extended TS List]
Ultrasound Image Storage	1.2.840.10008.5.1.4.1.1.6.1	[Extended TS List]
Enhanced US Volume Storage	1.2.840.10008.5.1.4.1.1.6.2	[Extended TS List]
Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7	[Extended TS List]
Multi-frame Single Bit Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7.1	[Extended TS List]
Multi-frame Grayscale Byte Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7.2	[Extended TS List]
Multi-frame Grayscale Word Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7.3	[Extended TS List]
Multi-frame True Color Secondary Capture Image Storage	1.2.840.10008.5.1.4.1.1.7.4	[Extended TS List]
Standalone Overlay Storage (Retired)	1.2.840.10008.5.1.4.1.1.8	[Default TS List]
Standalone Curve Storage (Retired)	1.2.840.10008.5.1.4.1.1.9	[Default TS List]
12-lead ECG Waveform Storage	1.2.840.10008.5.1.4.1.1.9.1.1	[Default TS List]
General ECG Waveform Storage	1.2.840.10008.5.1.4.1.1.9.1.2	[Default TS List]
Ambulatory ECG Waveform Storage	1.2.840.10008.5.1.4.1.1.9.1.3	[Default TS List]
Hemodynamic Waveform Storage	1.2.840.10008.5.1.4.1.1.9.2.1	[Default TS List]
Cardiac Electrophysiology Waveform Storage	1.2.840.10008.5.1.4.1.1.9.3.1	[Default TS List]
Basic Voice Audio Waveform Storage	1.2.840.10008.5.1.4.1.1.9.4.1	[Default TS List]

General Audit Waveform Storage	1.2.840.10008.5.1.4.1.1.9.4.2	[Default TS List]
Arterial Pulse Waveform Storage	1.2.840.10008.5.1.4.1.1.9.5.1	[Default TS List]
Respiratory Waveform Storage	1.2.840.10008.5.1.4.1.1.9.6.1	[Default TS List]
Standalone Modality LUT Storage (Retired)	1.2.840.10008.5.1.4.1.1.10	[Default TS List]
Standalone VOI LUT Storage (Retired)	1.2.840.10008.5.1.4.1.1.11	[Default TS List]
Grayscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.1	[Default TS List]
Color Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.2	[Default TS List]
Pseudo-Color Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.3	[Default TS List]
Blending Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.4	[Default TS List]
XA/XRF Grayscale Softcopy Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.5	[Default TS List]
Grayscale Planar MPR Volumetric Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.6	[Default TS List]
Compositing Planar MPR Volumetric Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.7	[Default TS List]
Advanced Blending Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.8	[Default TS List]
Volume Rendering Volumetric Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.9	[Default TS List]
Segmented Volume Rendering Volumetric Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.10	[Default TS List]
Multiple Volume Rendering Volumetric Presentation State Storage	1.2.840.10008.5.1.4.1.1.11.11	[Default TS List]
X-Ray Angiographic Image Storage	1.2.840.10008.5.1.4.1.1.12.1	[Extended TS List]
Enhanced XA Image Storage	1.2.840.10008.5.1.4.1.1.12.1.1	[Extended TS List]
X-Ray Radiofluoroscopic Image Storage	1.2.840.10008.5.1.4.1.1.12.2	[Extended TS List]
Enhanced XRF Image Storage	1.2.840.10008.5.1.4.1.1.12.2.1	[Extended TS List]
<i>X-Ray Angiographic Bi-Plane Image Storage (Retired)</i>	1.2.840.10008.5.1.4.1.1.12.3	[Extended TS List]
X-Ray 3D Angiographic Image Storage	1.2.840.10008.5.1.4.1.1.13.1.1	[Extended TS List]
X-Ray 3D Craniofacial Image Storage	1.2.840.10008.5.1.4.1.1.13.1.2	[Extended TS List]
Breast Tomosynthesis Image Storage	1.2.840.10008.5.1.4.1.1.13.1.3	[Extended TS List]
Breast Projection X-Ray Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.13.1.4	[Extended TS List]
Breast Projection X-Ray Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.13.1.5	[Extended TS List]
Intravascular Optical Coherence Tomography Image Storage – For Presentation	1.2.840.10008.5.1.4.1.1.14.1	[Extended TS List]

Intravascular Optical Coherence Tomography Image Storage – For Processing	1.2.840.10008.5.1.4.1.1.14.2	[Extended TS List]
Nuclear Medicine Image Storage	1.2.840.10008.5.1.4.1.1.20	[Extended TS List]
Parametric Map Storage	1.2.840.10008.5.1.4.1.1.30	[Extended TS List]
Raw Data Storage	1.2.840.10008.5.1.4.1.1.66	[Extended TS List]
Spatial Registration Storage	1.2.840.10008.5.1.4.1.1.66.1	[Extended TS List]
Spatial Fiducials Storage	1.2.840.10008.5.1.4.1.1.66.2	[Extended TS List]
Deformable Spatial Registration Storage	1.2.840.10008.5.1.4.1.1.66.3	[Extended TS List]
Segmentation Storage	1.2.840.10008.5.1.4.1.1.66.4	[Extended TS List]
Surface Segmentation Storage	1.2.840.10008.5.1.4.1.1.66.5	[Extended TS List]
Tractography Results Storage	1.2.840.10008.5.1.4.1.1.66.6	[Extended TS List]
Real World Value Mapping Storage	1.2.840.10008.5.1.4.1.1.67	[Extended TS List]
Surface Scan Mesh Storage	1.2.840.10008.5.1.4.1.1.68.1	[Extended TS List]
Surface Scan Point Cloud Storage	1.2.840.10008.5.1.4.1.1.68.2	[Extended TS List]
<i>VL Image Storage (Retired)</i>	1.2.840.10008.5.1.4.1.1.77.1	[Extended TS List]
VL Endoscopic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.1	[Extended TS List]
Video Endoscopic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.1.1	[Extended TS List]
VL Microscopic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.2	[Extended TS List]
Video Microscopic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.2.1	[Extended TS List]
VL Slide-Coordinates Microscopic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.3	[Extended TS List]
VL Photographic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.4	[Extended TS List]
Video Photographic Image Storage	1.2.840.10008.5.1.4.1.1.77.1.4.1	[Extended TS List]
Ophthalmic Photography 8 Bit Image Storage	1.2.840.10008.5.1.4.1.1.77.1.5.1	[Extended TS List]
Ophthalmic Photography 16 Bit Image Storage	1.2.840.10008.5.1.4.1.1.77.1.5.2	[Extended TS List]
Stereometric Relationship Storage	1.2.840.10008.5.1.4.1.1.77.1.5.3	[Extended TS List]
Ophthalmic Tomography Image Storage	1.2.840.10008.5.1.4.1.1.77.1.5.4	[Extended TS List]

Wide Field Ophthalmic Photography Stereographic Projection Image Storage	1.2.840.10008.5.1.4.1.1.77.1.5.5	[Extended TS List]
Wide Field Ophthalmic Photography 3D Coordinates Image Storage	1.2.840.10008.5.1.4.1.1.77.1.5.6	[Extended TS List]
Ophthalmic Optimal Coherence Tomography En Face Image Storage	1.2.840.10008.5.1.4.1.1.77.1.5.7	[Extended TS List]
Ophthalmic Optimal Coherence Tomography B-scan Volume Analysis Storage	1.2.840.10008.5.1.4.1.1.77.1.5.8	[Extended TS List]
VL Whole Slide Microscopy Image Storage	1.2.840.10008.5.1.4.1.1.77.1.6	[Extended TS List]
Lensometry Measurements Storage	1.2.840.10008.5.1.4.1.1.78.1	[Default TS List]
Autorefractive Measurements Storage	1.2.840.10008.5.1.4.1.1.78.2	[Default TS List]
Keratometry Measurements Storage	1.2.840.10008.5.1.4.1.1.78.3	[Default TS List]
Subjective Refraction Measurements Storage	1.2.840.10008.5.1.4.1.1.78.4	[Default TS List]
Visual Acuity Measurements Storage	1.2.840.10008.5.1.4.1.1.78.5	[Default TS List]
Spectacle Prescription Report Storage	1.2.840.10008.5.1.4.1.1.78.6	[Default TS List]
Ophthalmic Axial Measurements Storage	1.2.840.10008.5.1.4.1.1.78.7	[Default TS List]
Intraocular Lens Calculations Storage	1.2.840.10008.5.1.4.1.1.78.8	[Default TS List]
Macular Grid Thickness and Volume Report	1.2.840.10008.5.1.4.1.1.79.1	[Default TS List]
Ophthalmic Visual Field Static Perimetry Measurements Storage	1.2.840.10008.5.1.4.1.1.80.1	[Default TS List]
Ophthalmic Thickness Map Storage	1.2.840.10008.5.1.4.1.1.81.1	[Default TS List]
Corneal Topography Storage	1.2.840.10008.5.1.4.1.1.82.1	[Default TS List]
Basic Text SR	1.2.840.10008.5.1.4.1.1.88.11	[Default TS List]
Enhanced SR	1.2.840.10008.5.1.4.1.1.88.22	[Default TS List]
Comprehensive SR	1.2.840.10008.5.1.4.1.1.88.33	[Default TS List]
Comprehensive 3D SR	1.2.840.10008.5.1.4.1.1.88.34	[Default TS List]
Extensible SR Storage	1.2.840.10008.5.1.4.1.1.88.35	[Default TS List]
Procedure Log	1.2.840.10008.5.1.4.1.1.88.40	[Default TS List]
Mammography CAD SR	1.2.840.10008.5.1.4.1.1.88.50	[Default TS List]
Key Object Selection	1.2.840.10008.5.1.4.1.1.88.59	[Default TS List]
Chest CAD SR	1.2.840.10008.5.1.4.1.1.88.65	[Default TS List]
X-Ray Radiation Dose SR	1.2.840.10008.5.1.4.1.1.88.67	[Default TS List]

Radiopharmaceutical Radiation Dose SR	1.2.840.10008.5.1.4.1.1.88.68	[Default TS List]
Colon CAD SR	1.2.840.10008.5.1.4.1.1.88.69	[Default TS List]
Implantation Plan SR Document Storage	1.2.840.10008.5.1.4.1.1.88.70	[Default TS List]
Acquisition Context SR Storage	1.2.840.10008.5.1.4.1.1.88.71	[Default TS List]
Simplified Adult Echo SR Storage	1.2.840.10008.5.1.4.1.1.88.72	[Default TS List]
Patient Radiation Dose SR Storage	1.2.840.10008.5.1.4.1.1.88.73	[Default TS List]
Planning Imaging Agent Administration SR Storage	1.2.840.10008.5.1.4.1.1.88.74	[Default TS List]
Performed Imaging Agent Administration SR Storage	1.2.840.10008.5.1.4.1.1.88.75	[Default TS List]
Content Assessment Results Storage	1.2.840.10008.5.1.4.1.1.90.1	[Default TS List]
Encapsulated PDF Storage	1.2.840.10008.5.1.4.1.1.104.1	[Default TS List]
Encapsulated CDA Storage	1.2.840.10008.5.1.4.1.1.104.2	[Default TS List]
Encapsulated STL Storage	1.2.840.10008.5.1.4.1.1.104.3	[Default TS List]
Positron Emission Tomography Image Storage	1.2.840.10008.5.1.4.1.1.128	[Extended TS List]
Legacy Converted Enhanced PET Image Storage	1.2.840.10008.5.1.4.1.1.128.1	[Default TS List]
Standalone PET Curve Storage (Retired)	1.2.840.10008.5.1.4.1.1.129	[Default TS List]
Enhanced PET Image Storage	1.2.840.10008.5.1.4.1.1.130	[Extended TS List]
Basic Structured Display Storage	1.2.840.10008.5.1.4.1.1.131	[Extended TS List]
CT Performed Procedure Protocol Storage	1.2.840.10008.5.1.4.1.1.200.2	[Default TS List]
RT Image Storage	1.2.840.10008.5.1.4.1.1.481.1	[Extended TS List]
RT Dose Storage	1.2.840.10008.5.1.4.1.1.481.2	[Default TS List]
RT Structure Set Storage	1.2.840.10008.5.1.4.1.1.481.3	[Default TS List]
RT Beams Treatment Record Storage	1.2.840.10008.5.1.4.1.1.481.4	[Default TS List]
RT Plan Storage	1.2.840.10008.5.1.4.1.1.481.5	[Default TS List]
RT Brachy Treatment Record Storage	1.2.840.10008.5.1.4.1.1.481.6	[Default TS List]
RT Treatment Summary Record Storage	1.2.840.10008.5.1.4.1.1.481.7	[Default TS List]
RT Ion Plan Storage	1.2.840.10008.5.1.4.1.1.481.8	[Default TS List]
RT Ion Beams Treatment Record Storage	1.2.840.10008.5.1.4.1.1.481.9	[Default TS List]
RT Physician Intent Storage	1.2.840.10008.5.1.4.1.1.481.10	[Default TS List]
RT Segment Annotation Storage	1.2.840.10008.5.1.4.1.1.481.11	[Default TS List]
RT Beams Delivery Instruction Storage	1.2.840.10008.5.1.4.34.7	[Default TS List]

RT Brachy Application Setup Delivery Instruction Storage	1.2.840.10008.5.1.4.34.10	[Default TS List]
--	---------------------------	-------------------

2.2.1 Vitrea Connection Storage Server AE Specification

2.2.1.1 SOP Classes

Vitrea Connection Storage Server Application Entity provides Standard Conformance to the following SOP Classes:

SOP Classes for Vitrea Connection Storage Server AE

SOP Class Name	SOP Class UID	SCU	SCP
Verification SOP Class	1.2.840.10008.1.1	Yes	Yes
Patient Root Q/R Information Model - FIND	1.2.840.10008.5.1.4.1.2.1.1	Yes	No
Patient Root Q/R Information Model - MOVE	1.2.840.10008.5.1.4.1.2.1.2	Yes	No
Study Root Q/R Information Model - FIND	1.2.840.10008.5.1.4.1.2.2.1	Yes	No
Study Root Q/R Information Model - MOVE	1.2.840.10008.5.1.4.1.2.2.2	Yes	No
Storage Commitment Push Model	1.2.840.10008.5.1.20.1	No	Yes
[Default SOP Class List]		Yes	Yes

2.2.1.2 Association Establishment Policy

2.2.1.2.1 General

The Vitrea Connection Storage Server AE can both accept and propose Association Requests. The Vitrea Connection Storage Server AE will accept and propose Association Requests for the Verification, and Storage Services. It will accept Requests for Storage Commitment Push Model Services.

The DICOM standard application context name for DICOM 3.0 is always accepted and proposed:

DICOM application context name for Vitrea Connection Storage Server AE

Application Context Name	1.2.840.10008.3.1.1.1
--------------------------	-----------------------

2.2.1.2.2 Number of Associations

Vitrea Connection Storage Server can support multiple simultaneous Associations requested by peer AEs. Default is 10. This value can be configured through the configuration GUI.

Number of Associations accepted by Vitrea Connection Storage Server AE

Maximum number of simultaneous Associations	10 (Configurable)
---	-------------------

2.2.1.2.3 Asynchronous Nature

Vitrea Connection Storage Server does not support asynchronous communication. Multiple outstanding transactions are not supported. It allows up to one invoked and one performed operation on an Association (it is synchronous).

Asynchronous Nature as SCP for Vitrea Connection Storage Server AE

Maximum number of outstanding asynchronous transactions	1 (Not Configurable)
---	----------------------

2.2.1.2.4 Implementation Identifying Information

The implementation information for this Application Entity is:

DICOM Implementation Class and Version for Vitrea Connection Storage Server AE

Implementation Class UID	2.16.124.113638.2
Implementation Version Name	rialto-dcm4che2

2.2.1.3 Association Acceptance Policy

2.2.1.3.1 Activity - Receive Images

Description and Sequencing of Activities

A remote peer DICOM Application Entity, acting as a Storage SCU, establishes an association with Vitrea Connection Storage Server that accepts these Associations for the purpose of receiving supported SOP Class Instances requests.

In the default configuration any Calling and Called AET will be accepted. If the Called AET does not correspond to the actual Storage Server AET, only a Presentation Context for the Verification SOP Class will be accepted and the SCU can only verify the DICOM Association but cannot invoke any other related DICOM service.

The Storage Server AE may reject Association attempts as shown in the table below. The Result, Source and Reason / Diag columns represent the values returned in the corresponding fields of an ASSOCIATE-RJ PDU (see PS 3.8).

Association Rejection Reasons

Result	Source	Reason / Diag	Description
2 - rejected transient	provider	2 - local-limit exceeded	The (configurable) maximum number of simultaneous associations has been reached. An association request with the same parameters may succeed at a later time.
1 - rejected permanent	user	2 - application context name not supported	The Association request contained an unsupported Application Context Name. An association request with the same parameters will NOT succeed at a later time.
1 - rejected permanent	user	3 – calling AE title not recognized	The Association request contained an unrecognized Calling AE Title. An association request with the same parameters will NOT succeed at a later time unless configuration changes are made.

Accepted Presentation Contexts

Vitrea Connection Storage Server AE will accept Presentation Contexts as shown in the following table:

Accepted Presentation Contexts for Vitrea Connection Storage Server AE

Abstract Syntax		Transfer Syntax	Role	Ext. Neg.
Name	UID			

Verification SOP Class	1.2.840.10008.1.1	[Default TS List]	SCP	None
Storage Commitment Push Model	1.2.840.10008.1.20.1	[Default TS List]	SCP	None
[Default SOP Class List]			SCP	None

The preferred (default) Transfer Syntax is Explicit VR Little Endian.

SOP Specific Conformance

Specific Conformance for the Verification SOP Class

Vitreia Connection Storage Server provides standard conformance to the DICOM Verification Service Class as an SCP. The status code for the C-ECHO is described in the following table:

Vitreia Connection Storage Server C-Echo Response Status

Service Status	Further Meaning	Error Code	Reason
Success	Success	0000	The C-ECHO request is accepted

Specific Conformance for the Storage SOP Classes

The associated activity with the Storage service is the storage of DICOM data received over the network on a designated storage repository. The Vitrea Connection Storage Server AE will return a failure status if it is unable to store the received instance(s).

The Vitrea Connection Storage Server AE does not have any dependencies on the number of Associations used to send images to it. Images belonging to more than one Study or Series can be sent over a single or multiple Associations. Images belonging to a single Study or Series can also be sent over different Associations. There is no limit on either the number of SOP Instances or the maximum amount of total SOP Instance data that can be transferred over a single Association.

The Vitrea Connection Storage Server AE always retains the original DICOM data in DICOM Part 10 compliant file format, which is stored in the local storage system. The Vitrea Connection Storage Server AE is Level 2 (Full) conformant as a Storage SCP. In addition, all Private and SOP Class Extended Elements are maintained in the DICOM format files. In addition to saving all Elements in files, a subset of the Elements is stored in the Vitrea Connection Storage Server database to support query and retrieval requests.

By default, the Vitrea Connection Storage Server AE will not transcode information from one transfer syntax to another. All information is stored as it is received. However, it is possible to configure the Vitrea Connection Storage Server AE to effect transcoding between transfer syntaxes when storing images (as outlined in the table below). If the Vitrea Connection Storage Server AE is unable to transcode an image that is being stored, the Vitrea Connection Storage Server AE will store the image in its original format.

Transcode from Transfer Syntax	Transcode to Transfer Syntax
1.2.840.10008.1.2 (Implicit VR Little Endian)	1.2.840.10008.1.2.4.80 (JPEG-LS Lossless Image Compression)
1.2.840.10008.1.2.1 (Explicit VR Little Endian)	1.2.840.10008.1.2.4.80 (JPEG-LS Lossless Image Compression)
1.2.840.10008.1.2 .2 (Explicit VR Big Endian)	1.2.840.10008.1.2.4.80 (JPEG-LS Lossless Image Compression)

As an option, the Vitrea Connection Storage Server AE can be configured to update the metadata associated with the DICOM data based on information found in scheduled procedure steps (i.e. modality worklist reconciliation). The original DICOM data is not modified, but only the metadata. A complete history of metadata is stored, so that the updates performed can be examined and even rolled back by an administrator.

The Vitrea Connection Storage Server AE can be configured to selectively apply KOS rejection notes to subsequent DICOM C-Find and C-Move requests, depending on which AE Title is queried, as per the specification of Information Object Change Management.

Storage Server AE C-STORE Response

Vitrea Connection Storage Server C-STORE Response Status

Service Status	Further Meaning	Error Code	Reason
Success	Success	0000	The Composite SOP Instance was successfully received, verified, and stored in the system repository.
Error	Processing Failure	0110	This status is returned due to internal errors such as a processing failure response from the Vitrea Connection database or a file system operation. The appropriate Status will be sent in the C-STORE Response. Error indication message is output to the Service Log.
Error	Duplicate SOP Instance UID	D000	This status is returned if the instance already exists in the database and the Vitrea Connection Storage Server AE is configured to refuse duplicate instances. The appropriate Status will be sent in the C-STORE Response. Error indication message is output to the Service Log.
Warning	Data Set does not match SOP class	B007	This status is returned if the C-STORE Request specifies Attributes that are not specific as part of the Storage SOP class. Image transmission is considered successful. The appropriate SUCCESS Status will be sent in the C-STORE Response. Warning indication message is output to the Service Log.

Vitrea Connection Storage Server Service Communication Failure Behavior

Exception	Behavior
Timeout expiry for an expected DICOM Message Request (DIMSE level timeout). I.e., The Storage Server SCP AE is waiting for the next C-STORE Request on an open Association but the timer expires.	The Association is aborted by issuing a DICOM A-ABORT. Error message is output to the Service Log. If some Composite SOP Instances have already been successfully received, they are maintained in the database. They are not automatically discarded because of a later failure. Error message is output to the Service Log.
Timeout expiry for an expected DICOM PDU or TCP/ IP packet (Low-level timeout). I.e., The Storage Server SCP AE is waiting for the next C-STORE Data Set PDU but the timer expires.	The Association is aborted by issuing a DICOM A-ABORT. Error message is output to the Service Log. If a C-STORE Data Set has not been fully received, the data already received is discarded. If some Composite SOP Instances have already been successfully received over the Association, they are maintained in the database.

	Error message is output to the Service Log.
Association aborted by the SCU or the network layers indicate communication loss (i.e. low-level TCP/IP socket closure)	Error message is output to the Service Log. If some Composite SOP Instances have already been successfully received, they are maintained in the database. They are not automatically discarded because of a later failure.

Specific Conformance for the Storage Commitment SOP Class

The associated activity with the Storage Commitment Push Model service is the communication by the Vitrea Connection Storage Server AE to peer AEs that it has committed to permanently store Composite SOP Instances that have been sent to it. It thus allows peer AEs to determine whether the Vitrea Connection Storage Server has taken responsibility for the archiving of specific SOP Instances so that they can be flushed from the peer AE system.

Upon receipt of a Storage Commitment Push Model N-ACTION Request, the Vitrea Connection Storage Server AE provides an N-ACTION response/acknowledgement on the same association. The Vitrea Connection Storage Server AE takes the list of Composite SOP Instance UIDs specified in the N-ACTION Request and checks if they are present in the Vitrea Connection query/retrieve Server repository. As long as the Composite SOP Instance UIDs are present in the database, the Vitrea Connection Storage Server AE will consider those Composite SOP Instance UIDs to be successfully archived. The Vitrea Connection Storage Server AE does not require the Composite SOP Instances to actually be successfully written to archive media in order to commit to responsibility for maintaining these SOP Instances.

Once the Vitrea Connection Storage Server AE has checked for the existence of the specified Composite SOP Instances, it will then attempt to send the Notification request (N-EVENT-REPORT-RQ). By default, the Vitrea Connection Storage Server AE will request a new Association with the peer AE that made the original N-ACTION Request, but it may also be configured to provide a synchronous N-EVENT-REPORT on the initial association.

The Vitrea Connection Storage Server AE will not cache Storage Commitment Push Model N-ACTION Requests that specify Composite SOP Instances that have not yet been transferred to the Vitrea Connection Query/Retrieve Server. If a peer AE sends a Storage Commitment Push Model N-ACTION Request before the specified Composite SOP Instances are later sent over the same Association, the Vitrea Connection Storage Server AE will not commit to responsibility for such SOP Instances.

The Vitrea Connection Storage Server AE does not support the optional Storage Media File-Set ID & UID attributes in the N-ACTION request.

The Vitrea Connection Storage Server AE will support Storage Commitment Push Model requests for SOP Instances of any of the Storage SOP Classes that are also supported by the Vitrea Connection Storage Server AE:

Supported Referenced SOP Classes In Storage Commitment Push Model N-ACTION Requests

SOP Class Name	SOP Class UID
[Default SOP Class List]	

Storage Server AE N-ACTION Response

The Vitrea Connection Storage Server AE will return the following Status Code values in N-ACTION Responses:

Vitrea Connection Storage Server N-ACTION Response Status

Service Status	Further Meaning	Error Code	Reason
Success	Success	0000	The SCP has successfully received the Storage Commitment Push Model N-ACTION Request and can process the commitment request for the indicated SOP Instances.
Error	Processing Failure	0110	Indicates that the Storage Commitment Push Model N-

			ACTION Request cannot be parsed or fully processed due to a database or system failure.
Error	Missing Attribute	0120	Indicates that the Storage Commitment Push Model N-ACTION Request cannot be processed because a required attribute is missing from the N-ACTION Request Data Set.
Error	Missing Attribute Value	0121	Indicates that the Storage Commitment Push Model N-ACTION Request cannot be processed because a Type 1 attribute in the NACTION Request Data Set does not specify a value.

Storage Server AE N-EVENT Response

The Vitrea Connection Storage Server AE will exhibit the following Behavior according to the Status Code value returned in an N-EVENT-REPORT Response from a destination Storage Commitment Push Model SCU:

Vitrea Connection Storage Server N-EVENT Report Response Status

Service Status	Further Meaning	Error Code	Reason
Success	Success	0000	The SCU has successfully received the Storage Commitment Push Model N-EVENT-REPORT Request. Success indication message is output to the Service Log.
Warning	Attribute List Error	0107	Transmission of Storage Commitment Push Model N-EVENT-REPORT Request is considered successful. Warning indication message is output to the Service Log.
*	*	Any Other Code	This is treated as a permanent Failure. Error indication message is output to the Service Log.

The next table reports the Vitrea Connection Storage Server behavior in case of communication failure:

Vitrea Connection Storage Server AE Storage Commitment Push Model Communication Failure Behavior

Exception	Behavior
Timeout expiry for an expected DICOM Message Request (DIMSE level timeout). I.e., The Storage Server SCP AE is waiting for the next N-ACTION Request on an open Association but the timer expires.	The Association is aborted by issuing a DICOM A-ABORT. If some Composite SOP Instances have been successfully received over the same Association via the Storage Service, they are maintained in the database. They are not automatically discarded because of a later Storage Commitment messaging failure. Any previously received Storage Commitment Push Model N-ACTION Requests will still be fully processed. Error indication message is output to the Service Log.
Timeout expiry for an expected DICOM PDU or TCP/ IP packet (Low-level timeout).	The Association is aborted by issuing a DICOM A-ABORT. If some Composite SOP Instances have been successfully received over the same Association via the Storage Service, they are maintained in the database. They are not automatically discarded because of a later Storage Commitment messaging

	<p>failure.</p> <p>Any previously received Storage Commitment Push Model N-ACTION Request will still be fully processed.</p> <p>Error indication message is output to the Service Log.</p>
<p>Association A-ABORTed by the SCU or the network layers indicate communication loss (i.e., low-level TCP/IP socket closure).</p>	<p>The TCP/IP socket is closed.</p> <p>If some Composite SOP Instances have been successfully received over the same Association via the Storage Service, they are maintained in the database. They are not automatically discarded because of a later Storage Commitment messaging failure.</p> <p>Any previously received Storage Commitment Push Model N-ACTION Requests will still be fully processed.</p> <p>Error indication message is output to the Service Log.</p>

2.2.1.4 Association Initiation Policy

2.2.1.4.1 Activity – Synchronize Objects with an External Peer AE

Description and Sequencing of Activities

The Vitrea Connection Storage Server AE will initiate a new Association when attempting to synchronize objects that have been previously stored by an external Peer. The Storage Server AE will send the Association Request to the specific Peer, and upon successful negotiation of the required Presentation Context, a query request and possibly a retrieve request will be sent to the Peer.

Proposed Presentation Contexts

A list of Transfer Syntaxes can be configured for each AET. Vitrea Connection Storage Server AE will propose Presentation Contexts as shown in the following table:

Abstract Syntax		Transfer Syntax	Role	Ext. Neg.
Name	UID			
Verification SOP Class	1.2.840.10008.1.1	[Default TS List]	SCU	None
Patient Root Q/R Information Model - FIND	1.2.840.10008.5.1.4.1.2.1.1	[Default TS List]	SCU	None
Patient Root Q/R Information Model - MOVE	1.2.840.10008.5.1.4.1.2.1.2	[Default TS List]	SCU	None
Study Root Q/R Information Model - FIND	1.2.840.10008.5.1.4.1.2.2.1	[Default TS List]	SCU	None
Study Root Q/R Information Model - FIND	1.2.840.10008.5.1.4.1.2.2.2	[Default TS List]	SCU	None

SOP Specific Conformance

Specific Conformance for Verification SOP Class

Standard conformance is provided to the DICOM Verification Service Class as an SCU. The Verification Service as an SCU is only supported as a diagnostic service tool for network communication issues.

2.2.1.4.2 Activity – Archive Objects to an External Peer AE

Description and Sequencing of Activities

When configured as a cache, the Vitrea Connection Storage Server AE will initiate a new Association when attempting to archive objects that have been previously stored. The Storage Server AE will send the Association Request to the Peer that is serving as an archive, and upon successful negotiation of the required Presentation Context, a verification request, a storage request and possibly a storage commitment request will be sent to the Peer.

The Vitrea Connection Storage Server AE can be configured to export KOS rejection notes alongside archive objects when appropriately queried by external peer AEs.

Proposed Presentation Contexts

A list of Transfer Syntaxes can be configured for each AET. Vitrea Connection Storage Server AE will propose Presentation Contexts as shown in the following table:

Abstract Syntax		Transfer Syntax	Role	Ext. Neg.
Name	UID			
Verification SOP Class	1.2.840.10008.1.1	[Default TS List]	SCU	None
[Default SOP Class List]		[Default TS List]	SCU	None

2.2.2 Vitrea Connection Query/Retrieve Server AE Specification

2.2.2.1 SOP Classes

Vitrea Connection Query/Retrieve Server Application Entity provides Standard Conformance to the following SOP Classes:

SOP Classes for Vitrea Connection Query/Retrieve Server AE

SOP Class Name	SOP Class UID	SCU	SCP
Verification SOP Class	1.2.840.10008.1.1	Yes	Yes
Patient Root Q/R Information Model - FIND	1.2.840.10008.5.1.4.1.2.1.1	No	Yes
Patient Root Q/R Information Model - MOVE	1.2.840.10008.5.1.4.1.2.1.2	No	Yes
Study Root Q/R Information Model - FIND	1.2.840.10008.5.1.4.1.2.2.1	No	Yes
Study Root Q/R Information Model - MOVE	1.2.840.10008.5.1.4.1.2.2.2	No	Yes
[Default SOP Class List]		Yes	No

2.2.2.2 Association Establishment Policy

2.2.2.2.1 General

The Vitrea Connection Query/Retrieve Server AE can both accept and propose Association Requests. The Vitrea Connection Query/Retrieve Server AE will accept Association Requests for the Verification and Query/Retrieve Services. It will propose Associations for Verification and Storage Services.

The DICOM standard application context name for DICOM 3.0 is always accepted and proposed:

DICOM application context name for Vitrea Connection Query/Retrieve Server AE

Application Context Name	1.2.840.10008.3.1.1.1
--------------------------	-----------------------

2.2.2.2.2 Number of Associations

Vitrea Connection Query/Retrieve Server can support multiple simultaneous Associations requested by peer AEs. Default is 10. This value can be configured through the attribute "MaxClients" in the Application Server configuration.

Number of Associations accepted for Vitrea Connection Query/Retrieve Server AE

Maximum number of simultaneous Associations	10 (Configurable)
---	-------------------

2.2.2.2.3 Asynchronous Nature

Vitrea Connection Query/Retrieve Server does not support asynchronous communication. Multiple outstanding transactions are not supported. It allows up to one invoked and one performed operation on an Association (it is synchronous).

Asynchronous Nature as SCP for Vitrea Connection Query/Retrieve Server AE

Maximum number of outstanding asynchronous transactions	1 (Not Configurable)
---	----------------------

2.2.2.2.4 Implementation Identifying Information

The implementation information for this Application Entity is:

DICOM Implementation Class and Version for Vitrea Connection Query/Retrieve Server AE

Implementation Class UID	2.16.124.113638.2
Implementation Version Name	rialto-dcm4che2

2.2.2.3 Association Initiation Policy

2.2.2.3.1 Activity - Send Images Requested by an External Peer AE

Description and Sequencing of Activities

The Vitrea Connection Query/Retrieve AE will initiate a new Association when a valid C-MOVE Request has been performed by an external Peer. The Query/Retrieve AE will send the Association Request to the specific C-MOVE destination, and upon successful negotiation of the required Presentation Context the image transfer is started. In all cases an attempt will be made to transmit all the indicated images in a single Association but this may not always be possible. The Association will be released when all the images have been sent.

If an error occurs during transmission over an open Association, the image transfer is halted. The Query/Retrieve AE will not attempt to independently retry the image export.

By default, the Vitrea Connection Query/Retrieve Server AE will not transcode information from one transfer syntax to another when sending images to the specific C-MOVE destination. However, if the specific C-MOVE destination does not support the transfer syntax that the image was stored in, either as it was received originally or as transcoded by the Vitrea Connection Storage Server AE, the Vitrea Connection Query/Retrieve Server AE will attempt the following transcodings:

Transcode from Transfer Syntax	Transcode to Transfer Syntax
1.2.840.10008.1.2.4.80 (JPEG-LS Lossless Image Compression)	1.2.840.10008.1.2 (Implicit VR Little Endian)

Proposed Presentation Contexts

A list of Transfer Syntaxes can be configured for each AET. The default Transfer Syntaxes are the Transfer Syntax of the originally received instance plus the [Default TS List].

Vitrea Connection Query/Retrieve AE will propose Presentation Contexts as shown in the following table:

Proposed Presentation Contexts by the Vitrea Connection Query/Retrieve Server AE

Abstract Syntax			
-----------------	--	--	--

Name	UID	Transfer Syntax	Role	Ext. Neg.
Verification SOP Class	1.2.840.10008.1.1	[Default TS List]	SCU	None
[Default SOP Class List]			SCU	None

SOP Specific Conformance

Specific Conformance for Verification SOP Class

Standard conformance is provided to the DICOM Verification Service Class as an SCU. The Verification Service as an SCU is actually only supported as a diagnostic service tool for network communication issues.

Specific Conformance for Image, SR and KIN SOP Classes

The Query/Retrieve Server AE will exhibit the following behavior according to the Status Code value returned in a C-STORE Response from a destination C-STORE SCP:

Vitreia Connection Query/Retrieve Server AE C-STORE Response Status Handling Behavior

Service Status	Further Meaning	Error Code	Behavior
Success	Success	0000	The remote Storage SCP has successfully stored the exported SOP Instance. The appropriate PENDING or SUCCESS Status will be sent in the C-MOVE Response. Success indication message is output to the Service Log.
Refused	Out of Resources	A700 – A7FF	This is treated as a permanent Failure. The appropriate Status will be sent in the C-MOVE Response. Error indication message is output to the Service Log.
Error	Cannot Understand	C000 – CFFF	This is treated as a permanent Failure. The appropriate Status will be sent in the C-MOVE Response. Error indication message is output to the Service Log.
Warning	Coercion of Data Elements	B000	Image transmission is considered successful. The appropriate PENDING or SUCCESS Status will be sent in the C-MOVE Response. Success indication message is output to the Service Log.
Warning	Data Set does not match SOP Class	B007	Image transmission is considered successful. The appropriate PENDING or SUCCESS Status will be sent in the C-MOVE Response. Warning indication message is output to the Service Logs.
Warning	Elements Discarded	B006	Image transmission is considered successful. The appropriate PENDING or SUCCESS Status will be sent in the C-MOVE Response. Success indication message is output to the Service Log.
Warning	Attribute List Error	0107	Image transmission is considered successful.

			The appropriate PENDING or SUCCESS Status will be sent in the C-MOVE Response. Success indication message is output to the Service Log.
Warning	Attribute Value Out of Range	0116	Image transmission is considered successful. The appropriate PENDING or SUCCESS Status will be sent in the C-MOVE Response. Success indication message is output to the Service Log.
*	*	Any other status code	This is treated as a permanent Failure. The appropriate Status will be sent in the C-MOVE Response. Error indication message is output to the Service Log.

Vitrea Connection Query/Retrieve Server AE Communication Failure Behavior

Exception	Behavior
Timeout expiry for an expected DICOM Message Request (DIMSE level timeout).	The Association is aborted by issuing a DICOM A-ABORT. The appropriate Status will be sent in the CMOVE Response. Error indication message is output to the Service Logs.
Timeout expiry for an expected DICOM PDU or TCP/ IP packet (Low-level timeout).	The Association is aborted by issuing a DICOM A-ABORT. The appropriate Status will be sent in the CMOVE Response. Error indication message is output to the Service Logs.
Association A-ABORTed by the SCU or the network layers indicate communication loss (i.e. low-level TCP/IP socket closure).	The Association is aborted by issuing a DICOM A-ABORT. The appropriate Status will be sent in the CMOVE Response. Error indication message is output to the Service Logs.

2.2.2.3.2 Activity - Retrieve Images Requested by an External Peer AE

Description and Sequencing of Activities

The Vitrea Connection Query/Retrieve AE will initiate a new Association when a valid C-MOVE Request has been performed by an external Peer but one or more of the requested objects is not found in local storage. The Query/Retrieve AE will send an Association Request to the AE that originally stored the objects, and upon successful negotiation of the required Presentation Context a retrieve request is started. The Association will be released when all the images have been sent.

Proposed Presentation Contexts

A list of Transfer Syntaxes can be configured for each AET. The default Transfer Syntaxes are the Transfer Syntax of the originally received instance plus the [Default TS List].

The following Transfer Syntaxes are represented by the term [Default TS List]:

[Default TS List] Transfer Syntax List for Vitrea Connection Retrieve Server AE

Name List	UID List
Implicit VR Little Endian	1.2.840.10008.1.2
Explicit VR Little Endian	1.2.840.10008.1.2.1

Vitrea Connection Query/Retrieve AE will propose Presentation Contexts as shown in the following table:

Proposed Presentation Contexts by the Vitrea Connection Query/Retrieve Server AE

Abstract Syntax		Transfer Syntax	Role	Ext. Neg.
Name	UID			
Verification SOP Class	1.2.840.10008.1.1	[Default TS List]	SCU	None
Patient Root Q/R Information Model – MOVE	1.2.840.10008.5.1.4.1.2.1.2	[Default TS List]	SCU	Relational
Study Root Q/R Information Model - MOVE	1.2.840.10008.5.1.4.1.2.2.2	[Default TS List]	SCU	Relational

SOP Specific Conformance

Specific Conformance for Verification SOP Class

Standard conformance is provided to the DICOM Verification Service Class as an SCU. The Verification Service as an SCU is actually only supported as a diagnostic service tool for network communication issues.

2.2.2.4 Association Acceptance Policy

2.2.2.4.1 Activity - Handling Query and Retrieval Requests

Description and Sequencing of Activities

The Query/Retrieve Server AE accepts Associations only if they have valid Presentation Contexts. If none of the requested Presentation Contexts are accepted, the Association Request itself is rejected. It can be configured to only accept Associations from certain Application Entities.

If Query/Retrieve Server AE receives a query (C-FIND) request, the response(s) will be sent over the same Association used to send the C-FIND-Request.

If Query/Retrieve Server AE receives a retrieval (C-MOVE) request, the responses will be sent over the same Association used to send the C-MOVE-Request.

The Query/Retrieve Server AE will send the requested SOP Instances to the C-MOVE Destination over a newly created Association and report in the C-MOVE-Response any success or failure status of each attempt to send a Composite SOP Instance. The Query/Retrieve Server AE will not transcode information the stored transfer syntax to a different transfer syntax when sending the requested SOP Instances.

The Query/Retrieve Server AE may reject Association attempts as shown in the table below. The Result, Source and Reason / Diag columns represent the values returned in the corresponding fields of an ASSOCIATE-RJ PDU (see PS 3.8).

Association Rejection Reasons

Result	Source	Reason / Diag	Description
2 - rejected transient	provider	2 - local-limit exceeded	The (configurable) maximum number of simultaneous associations has been reached. An association request with the same parameters may succeed at a later time.
1 - rejected permanent	user	2 - application context name not supported	The Association request contained an unsupported Application Context Name. An association request with the same parameters will NOT succeed at a later time.
1 - rejected permanent	user	3 – calling AE title not recognized	The Association request contained an unrecognized Calling AE Title. An association request with the same parameters will NOT succeed at a later time unless configuration changes are made.

Accepted Presentation Contexts

Vitreia Connection Query/Retrieve Server AE will accept Presentation Contexts as shown in the following table:

Accepted Presentation Contexts for Vitreia Connection Query/Retrieve Server AE

Abstract Syntax		Transfer Syntax		Role	Ext. Neg.
Name	UID	Name List	UID List		
Verification SOP Class	1.2.840.10008.1.1	Implicit VR Little Endian Explicit VR Little Endian	1.2.840.10008.1.2 1.2.840.10008.1.2.1	SCP	None
Patient Root Q/R Information Model – FIND	1.2.840.10008.5.1.4.1.2.1.1	Implicit VR Little Endian Explicit VR Little Endian	1.2.840.10008.1.2 1.2.840.10008.1.2.1	SCP	None
Patient Root Q/R Information Model – MOVE	1.2.840.10008.5.1.4.1.2.1.2	Implicit VR Little Endian Explicit VR Little Endian	1.2.840.10008.1.2 1.2.840.10008.1.2.1	SCP	Relational
Study Root Q/R Information Model – FIND	1.2.840.10008.5.1.4.1.2.2.1	Implicit VR Little Endian Explicit VR Little Endian	1.2.840.10008.1.2 1.2.840.10008.1.2.1	SCP	None
Study Root Q/R Information Model – MOVE	1.2.840.10008.5.1.4.1.2.2.2	Implicit VR Little Endian Explicit VR Little Endian	1.2.840.10008.1.2 1.2.840.10008.1.2.1	SCP	Relational

If the Called AET is not corresponding to the actual Query/Retrieve Server AET, only the Presentation Context for the Verification SOP Class will be accepted.

SOP Specific Conformance

Specific Conformance for Verification SOP Class

Vitreia Connection Query/Retrieve Server provides standard conformance to the DICOM Verification Service Class as an SCP. The status code for the C-ECHO is described in the following table:

Vitreia Connection Query/Retrieve Server C-Echo Response Status

Service Status	Further Meaning	Error Code	Reason
Success	Success	0000	The C-ECHO request is accepted

Specific Conformance for Query SOP Classes

The Query/Retrieve Server AE supports hierarchical queries and not relational queries. There are no attributes, which are always returned by default. Only those attributes requested in the query identifier are returned.

Query responses always return values from the Query/Retrieve Server database. Exported SOP Instances are always updated with the latest values kept in the database prior to export. Thus, a change in Patient demographic information will be contained in both the C-FIND Responses and any Composite SOP Instances exported to a C-MOVE Destination AE.

If no Issuer of Patient ID and Issuer of Patient ID Qualifiers Sequence is present, default configurable values are assumed.

Patient Root Information Model

All required search keys on the Patient and Study levels are supported.

Study Root Information Model

All the required search keys on the Study level are supported.

Patient Root C-FIND SCP Supported Elements

SOP Common - Attribute Name	Tag	VR	Types of Matching
Specific Character Set	0008,0005	CS	NONE
Instance Availability	0008,0056	CS	S,U
Patient Level - Attribute Name	Tag	VR	Types of Matching
Patient's Name	0010,0010	PN	S,W,U
Patient ID	0010,0020	LO	S,W,U
Issuer of Patient ID	0010,0021	LO	S,U
Issuer of Patient ID Qualifiers Sequence	0010,0024	SQ	S,U
>Universal Entity ID	0040,0032	UT	S
>University Entity ID Type	0040,0033	CS	S
Other Patient IDs Sequence	0010,1002	SQ	U
Patient's Birth Date	0010,0030	DA	S,U,R
Patient's Birth Time	0010,0032	TM	S,U,R
Patient's Sex	0010,0040	CS	S,U
Study Level - Attribute Name	Tag	VR	Types of Matching
Study Instance UID	0020,000D	UI	S,L,U
Study ID	0020,0010	SH	S,W,U
Study Date	0008,0020	DA	S,W,U,R
Study Time	0008,0030	TM	S,W,U,R
Accession Number	0008,0050	SH	S,W,U
Modalities In Study	0008,0061	CS	U
Study Description	0008,1030	LO	U
Referring Physician's Name	0008,0090	PN	U
Institution Name	0008,0080	LO	S,W,U
Procedure Code Sequence	0008,1032	SQ	S,U

Number of Study Related Series	0020,1206	IS	U
Number of Study Related Instances	0020,1208	IS	U
SOP Classes in Study	0008,0062	UI	U
Series Level - Attribute Name	Tag	VR	Types of Matching
Series Instance UID	0020,000E	UI	S,U
Modality	0008,0060	CS	S,U
Series Date	0008,0021	DA	S,R,U
Series Time	0008,0031	TM	S,R,U
Series Description	0008,103E	LO	S,W,U
Laterality	0020,0060	CS	S,U
Anatomic Region Sequence	0008,2218	SQ	U
Body Part Examined	0018,0015	CS	S,W,U
Manufacturer	0008,0070	LO	S,W,U
Manufacturer's Model Name	0008,1090	LO	S,W,U
Frame of Reference UID	0020,0052	UI	S,U
Performed Procedure Step Description	0040,0254	LO	S,W,U
Protocol Name	0018,1030	LO	S,W,U
Station Name	0008,1010	SH	S,W,U
Institutional Department Name	0008,1040	LO	S,W,U
Performing Physician's Name	0008,1050	PN	S,W,U
Performed Procedure Step Start Date	0040,0244	DA	S,R,U
Performed Procedure Step Start Time	0040,0245	TM	S,R,U
Corrected Image	0028,0051	LO	S,W,U
Date of Last Calibration	0018,1200	DA	U
Device Serial Number	0018,1000	LO	S,U
Institution Address	0008,0081	ST	S,U
Institutional Department Name	0008,1040	LO	S,W,U
Institution Name	0008,0080	LO	S,W,U
Operators' Name	0008,1070	PN	S,W,U
Performed Procedure Step Status	0040,0252	CS	S,W,U

Pixel Padding Value	0028,0120	US or SS	S,U
Presentation Intent Type	0008,0068	CS	S,W,U
Series Number	0020,0011	IS	S,U
Series Type	0054,1000	CS	S,W,U
Software Version(s)	0018,1020	LO	S,W,U
Spatial Resolution	0018,1050	DS	S,W,U
Station Name	0008,1010	SH	S,W,U
Date of Last Calibration	0018,1200	DA	S,R,U
Time of Last Calibration	0018,1201	TM	S,R,U
Number of Series Related Instances	0020,1209	IS	U
SOP Level - Attribute Name	Tag	VR	Types of Matching
SOP Instance UID	0008,0018	UI	S,L,U
SOP Class UID	0008,0016	UI	S,U
Rows	0028,0010	US	S,U
Columns	0028,0011	US	S,U
Number of Frames	0028,0008	IS	S,U
Bits Allocated	0028,0100	US	S,U
Content Date	0008,0023	DA	S,R,U
Content Time	0008,0033	TM	S,R,U
Observation Date Time	0040,A032	DT	S,U
Concept Name Code Sequence	0040,A043	SQ	S,U
Content Label	0070,0080	CS	S,W,U
Content Description	0070,0081	LO	S,W,U
Presentation Creation Date	0070,0082	DA	S,R,U
Presentation Creation Time	0070,0083	TM	S,R,U
Content Creator's Name	0070,0084	PN	S,W,U
Transfer Syntax UID	0002,0010	UI	S,U
Requested Procedure ID	0040,1001	SH	S,U

Study Root C-FIND SCP Supported Elements

SOP Common - Attribute Name	Tag	VR	Types of Matching
[All attributes from SOP Common Level of Patient Root Model]			
Study Level - Attribute Name	Tag	VR	Types of Matching
[All attributes from Patient and Study Levels of Patient Root Model]			
Series Level - Attribute Name	Tag	VR	Types of Matching
[All attributes from Series Level of Patient Root Model]			
SOP Level - Attribute Name	Tag	VR	Types of Matching
[All attributes from SOP Level of Patient Root Model]			

Types of Matching:

The types of Matching supported by the C-FIND SCP. An "S" indicates the identifier attribute uses Single Value Matching, an "R" indicates Range Matching, a "W" indicates wildcard matching, a 'U' indicates Universal Matching (i.e. return keys), and an 'L' indicates that UID lists are supported for matching. NONE" means that single, range and wildcard matches are not supported on the key in question.

Cross-Enterprise Queries

If the C-FIND request includes an Issuer of Patient ID attribute, and Vitrea Connection is configured to recognize that Issuer as a Cross-Enterprise Issuer, Vitrea Connection will expand the query across all Issuers that the patient is known by (i.e. a cross-enterprise query) and will return the results of that expanded query.

Vitrea Connection Query/Retrieve Server C-FIND Response Status

Service Status	Further Meaning	Error Code	Reason
Success	Success	0000	Matching is complete. No final identifier is supplied.
Error	Processing Failure	0110	This status is returned due to internal errors such as a processing failure response on a query of the Vitrea Connection database. The appropriate Status will be sent in the C-FIND Response. Error indication message is output to the Service Log.
Cancel	Matching terminated due to Cancel Request	FE00	The C-FIND SCU sent a Cancel Request. This has been acknowledged and the search for matches has been halted.. Cancel indication message is output to the Service Log.
Pending	Matches are continuing and current match is supplied	FF00	Indicates that the search for further matches is continuing. This is returned when each successful match is returned and when further matches are forthcoming. Pending indication message is output to the Service Log.

Specific Conformance for Retrieve SOP Classes

The Query/Retrieve Server AE supports hierarchical and relational retrieves.

Exported SOP Instances are always updated with the latest values kept in the database prior to export.

Vitrea Connection Query/Retrieve Server C-MOVE Response Status

Service Status	Further Meaning	Error Code	Reason
Success	Success	0000	Matching is complete. No final identifier is supplied.
Error	Processing Failure	0110	<p>This status is returned due to internal errors such as a processing failure response on a query of the Vitrea Connection database.</p> <p>The appropriate Status will be sent in the C-MOVE Response.</p> <p>Error indication message is output to the Service Log.</p>
Error	Unable to perform sub-operations	A702	<p>C-STORE sub-operations cannot be performed due to failure of an Association Request or a C-STORE Request.</p> <p>Error indication message is output to the Service Log.</p>
Error	Move Destination unknown	A801	<p>The Destination Application Entity named in the C-MOVE Request is unknown to Query/Retrieve SCP AE.</p> <p>Error indication message is output to the Service Log.</p>
Error	Identifier does not match SOP Class	A900	<p>The C-MOVE identifier contains invalid Elements or values, or is missing mandatory Elements or values for the specified SOP Class or retrieval level.</p> <p>Error indication message is output to the Service Log.</p>
Error	Unable to process	Cxxx	<p>The Move Destination AET is missing in the C-MOVE Request.</p> <p>Error indication message is output to the Service Log.</p>
Pending	Sub-operations are continuing	FF00	<p>A Response with this Status Code is sent every time a Composite SOP Instance has been successfully sent to the C-MOVE Destination AE.</p>

Vitrea Connection Query/Retrieve Server Service Communication Failure Behavior

Exception	Behavior
Timeout expiry for an expected DICOM Message Request (DIMSE level timeout). I.e. the Query/Retrieve Server SCP AE is waiting for the next C-FIND or C-MOVE Request on an open Association but the timer expires.	The Association is aborted by issuing a DICOM A-ABORT. Error message is output to the Service Log.
Timeout expiry for an expected DICOM PDU or TCP/ IP packet (Low-level timeout). I.e. The Query/Retrieve Server AE is waiting for the next message PDU but the timer expires.	The Association is aborted by issuing a DICOM A-ABORT. Error message is output to the Service Log.
Association aborted by the SCU or the network layers indicate communication loss (i.e. low-level TCP/IP socket closure)	Error message is output to the Service Log.

2.2.3 Vitrea Connection Modality Worklist AE Specification

2.2.3.1 SOP Classes

Vitrea Connection MWL Server Application Entity provides Standard Conformance to the following SOP Classes:

SOP Classes for Vitrea Connection Modality Worklist AE

SOP Class Name	SOP Class UID	SCU	SCP
Verification SOP Class	1.2.840.10008.1.1	No	Yes
Modality Worklist Information Model - FIND	1.2.840.10008.5.1.4.31	No	Yes

2.2.3.2 Association Establishment Policy

2.2.3.2.1 General

The Vitrea Connection Modality Worklist AE can both accept and propose Association Requests. The Vitrea Connection Modality Worklist AE will accept Association Requests for the Verification and Modality Worklist Query Services.

The DICOM standard application context name for DICOM 3.0 is always accepted and proposed:

DICOM application context name for Vitrea Connection Query/Retrieve Server AE

Application Context Name	1.2.840.10008.3.1.1.1
--------------------------	-----------------------

2.2.3.2.2 Number of Associations

Vitrea Connection MWL Server can support multiple simultaneous Associations requested by peer AEs. Default is 10. This value can be configured through the attribute "MaxClients" in the Application Server configuration.

Number of Associations accepted for Vitrea Connection Modality Worklist AE

Maximum number of simultaneous Associations	10 (Configurable)
---	-------------------

2.2.3.2.3 Asynchronous Nature

Vitrea Connection MWL Server does not support asynchronous communication. Multiple outstanding transactions are not supported. It allows up to one invoked and one performed operation on an Association (it is synchronous).

Asynchronous Nature as SCP for Vitrea Connection Modality Worklist AE

Maximum number of outstanding asynchronous transactions	1 (Not Configurable)
---	----------------------

2.2.3.2.4 Implementation Identifying Information

The implementation information for this Application Entity is:

DICOM Implementation Class and Version for Vitrea Connection Modality Worklist AE

Implementation Class UID	2.16.124.113638.2
Implementation Version Name	rialto-dcm4che2

2.2.3.3 Association Acceptance Policy

2.2.3.3.1 Activity - Handling Query and Retrieval Requests

Description and Sequencing of Activities

The Modality Worklist AE accepts Associations only if they have valid Presentation Contexts. If none of the requested Presentation Contexts are accepted, the Association Request itself is rejected. It can be configured to only accept Associations from certain Application Entities.

If Modality Worklist AE receives a query (C-FIND) request, the response(s) will be sent over the same Association used to send the C-FIND-Request.

The Modality Worklist AE may reject Association attempts as shown in the table below. The Result, Source and Reason / Diag columns represent the values returned in the corresponding fields of an ASSOCIATE-RJ PDU (see PS 3.8).

Association Rejection Reasons

Result	Source	Reason / Diag	Description
2 - rejected transient	provider	2 - local-limit exceeded	The (configurable) maximum number of simultaneous associations has been reached. An association request with the same parameters may succeed at a later time.
1 - rejected permanent	user	2 - application context name not supported	The Association request contained an unsupported Application Context Name. An association request with the same parameters will NOT succeed at a later time.
1 - rejected permanent	user	3 – calling AE title not recognized	The Association request contained an unrecognized Calling AE Title. An association request with the same parameters will NOT succeed at a later time unless configuration changes are made.

Accepted Presentation Contexts

Vitrea Connection Query/Retrieve Server AE will accept Presentation Contexts as shown in the following table:

Accepted Presentation Contexts for Vitrea Connection Query/Retrieve Server AE

Abstract Syntax		Transfer Syntax		Role	Ext. Neg.
Name	UID	Name List	UID List		
Verification SOP Class	1.2.840.10008.1.1	Implicit VR Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		
Modality Worklist Information Model – FIND	1.2.840.10008.5.1.4.31	Implicit VR Little Endian	1.2.840.10008.1.2	SCP	None
		Explicit VR Little Endian	1.2.840.10008.1.2.1		

If the Called AET is not corresponding to the actual Modality Worklist AET, only the Presentation Context for the Verification SOP Class will be accepted.

SOP Specific Conformance

Specific Conformance for Verification SOP Class

Vitreia Connection MWL Server provides standard conformance to the DICOM Verification Service Class as an SCP. The status code for the C-ECHO is described in the following table:

Vitreia Connection Query/Retrieve Server C-Echo Response Status

Service Status	Further Meaning	Error Code	Reason
Success	Success	0000	The C-ECHO request is accepted

Specific Conformance for Query SOP Classes

The Modality Worklist AE supports modality worklist queries.

Query responses always return values from the MWL Server database

Modality Worklist C-FIND SCP Supported Elements

SOP Common - Attribute Name	Tag	VR	Types of Matching
Specific Character Set	0008,0005	CS	NONE
Scheduled Procedure Step - Attribute Name	Tag	VR	Types of Matching
Scheduled Procedure Step Sequence	0040,0100	SQ	U
>Scheduled Station AE Title	0040,0001	AE	S,U
>Scheduled Procedure Step Start Date	0040,0002	DA	S,R,U
>Scheduled Procedure Step Start Time	0040,0003	TM	S,R,U
>Modality	0008,0060	CS	S,U
>Scheduled Performing Physician's Name	0040,0007	PN	S,W,U
>Scheduled Station Name	0040,0010	SH	S,U
>Scheduled Procedure Step Location	0040,0011	SH	S,U
>Scheduled Protocol Code Sequence	0040,0008	SQ	U
>>Code Value	0008,0100	SH	U
>>Coding Scheme Designator	0008,0102	SH	U
>>Code Meaning	0008,0104	LO	U
>Pre-Medication	0040,0012	LO	U
>Scheduled Procedure Step ID	0040,0009	SH	U
>Requested Contrast Agent	0032,1070	LO	U
>Scheduled Procedure Step Status	0040,0020	CS	S,U
Requested Procedure - Attribute Name	Tag	VR	Types of Matching
Requested Procedure ID	0040,1001	SH	S,U

Reason for the Requested Procedure	0040,1002	LO	U
Reason for the Requested Procedure Code Sequence	0040,100A	SQ	U
>Code Value	0008,0100	SH	U
>Coding Scheme Designator	0008,0102	SH	U
>Code Meaning	0008,0104	LO	U
Requested Procedure Description	0032,1060	LO	U
Requested Procedure Code Sequence	0032,1064	SQ	U
>Code Value	0008,0100	SH	U
>Coding Scheme Designator	0008,0102	SH	U
>Code Meaning	0008,0104	LO	U
Study Instance UID	0020,000D	UI	S,U
Requested Procedure Priority	0040,1003	SH	U
Patient Transport Arrangements	0040,1004	LO	U
Imaging Service Request - Attribute Name	Tag	VR	Types of Matching
Accession Number	0008,0050	SH	S,U
Requesting Physician	0032,1032	PN	U
Referring Physician's Name	0008,0090	PN	U
Placer Order Number	0040,2016	LO	U
Filler Order Number	0040,2017	LO	U
Reason for Imaging Service Request	0040,2001	LO	U
Order Entered By	0040,2008	PN	U
Order Enterer's Location	0040,2009	SH	U
Order Callback Phone Number	0040,2010	SH	U
Visit Identification - Attribute Name	Tag	VR	Types of Matching
Admission ID	0038,0010	LO	U
Issuer of Admission ID	0038,0011	LO	U
Visit Status - Attribute Name	Tag	VR	Types of Matching
Current Patient Location	0030,0300	LO	U
Patient Identification - Attribute Name	Tag	VR	Types of Matching
Patient's Name	0010,0010	PN	S,W,U

Patient ID	0010,0020	LO	S,U
Issuer of Patient ID	0010,0021	LO	S,U
Ethnic Group	0010,2016	SH	U
Patient Demographics - Attribute Name	Tag	VR	Types of Matching
Patient's Birth Date	0010,0030	DA	S,R,U
Patient's Sex	0010,0040	CS	U
Confidentiality Constraint on Patient Data Description	0040,3001	LO	U
Region of Residence	0010,2152	LO	U
Military Rank	0010,1080	LO	U
Patient Medical - Attribute Name	Tag	VR	Types of Matching
Patient State	0038,0500	LO	U
Pregnancy Status	0010,21C0	US	U
Medical Alerts	0010,2000	LO	U
Allergies	0010,2110	LO	U

Types of Matching:

The types of Matching supported by the C-FIND SCP. An "S" indicates the identifier attribute uses Single Value Matching, an "R" indicates Range Matching, a "W" indicates wildcard matching, a 'U' indicates Universal Matching (i.e. return keys), and an 'L' indicates that UID lists are supported for matching. "NONE" means that single, range and wildcard matches are not supported on the key in question.

Vitrea Connection MWL Server C-FIND Response Status

Service Status	Further Meaning	Error Code	Reason
Success	Success	0000	Matching is complete. No final identifier is supplied.
Error	Processing Failure	0110	This status is returned due to internal errors such as a processing failure response on a query of the Vitrea Connection database. The appropriate Status will be sent in the C-FIND Response. Error indication message is output to the Service Log.
Cancel	Matching terminated due to Cancel Request	FE00	The C-FIND SCU sent a Cancel Request. This has been acknowledged and the search for matches has been halted. Cancel indication message is output to the Service Log.
Pending	Matches are continuing and current match is supplied	FF00	Indicates that the search for further matches is continuing. This is returned when each successful match is returned and when further matches are forthcoming. Pending indication message is output to the Service Log.

Vitreia Connection MWL Server Service Communication Failure Behavior

Exception	Behavior
Timeout expiry for an expected DICOM Message Request (DIMSE level timeout). I.e. the Query/Retrieve Server SCP AE is waiting for the next C-FIND or C-MOVE Request on an open Association but the timer expires.	The Association is aborted by issuing a DICOM A-ABORT. Error message is output to the Service Log.
Timeout expiry for an expected DICOM PDU or TCP/ IP packet (Low-level timeout). I.e. The Query/Retrieve Server AE is waiting for the next message PDU but the timer expires.	The Association is aborted by issuing a DICOM A-ABORT. Error message is output to the Service Log.
Association aborted by the SCU or the network layers indicate communication loss (i.e. low-level TCP/IP socket closure)	Error message is output to the Service Log.

2.2.4 Vitreia Connection STOW-RS SCP Specification

Vitreia Connection includes a STOW-RS SCP Application Entity that can be used to facilitate storage of DICOM objects via a web based interface. This service may be configured to be solely accessible using TLS if required. The following table describes the configurable attributes of the STOW-RS SCP.

STOW Configuration Parameters

Parameter	Configurable (Yes/No)	Default Value
Listening URL	Yes	http://localhost:8080/studies/{study}

The following table describes the supported MIME types of the STOW-RS Resource.

STOW MIME Types

MIME Type	Description
application/dicom	The incoming multipart message body content will be a "Part 10 File" that includes the meta-header.

STOW-RS Query Parameters

The STOW-RS service has no query parameters.

2.2.5 Vitreia Connection QIDO-RS Provider Specification

2.2.5.1 QIDO-RS Search for Studies

QIDO-RS Search for Studies

Parameter	Restrictions
Media Types	Restricted to 'multipart/related; type=application/dicom+xml' or 'application/json'
Matching Attributes	Refer QIDO-RS Study Attribute Matching
Limit and Offset supported	Yes
Fuzzy Person Name Matching supported	Yes

QIDO-RS Study Attribute Matching

Study Level - Attribute Name	Tag	VR	Types of Matching
Patient ID	0010.0020	LO	U,S,W
Issuer Of Patient ID	0010.0021	LO	U
Issuer Of Patient ID Qualifiers Sequence	0010.0024	SQ	U
Patient Name	0010.0010	PN	U,S,W
Patient Sex	0010.0040	CS	U
Patient Birth Date	0010.0030	DA	U
Study Instance UID	0020,000D	UI	U,S,L
Accession Number	0008,0050	SH	U,S,W
Issuer Of Accession Number Sequence	0008.0051	SQ	U
Study ID	0020,0010	SH	U,S,W
Study Date	0008,0020	DA	U,S,R
Study Time	0008,0030	TM	U,S,R
Modalities In Study	0008,0061	CS	U,S,W
Referring Physician's Name	0008,0090	PN	U,S,W
Study Description	0008,1030	LO	U
Instance Availability	0008.0056	CS	U
Institution Name	0008,0080	LO	U
Study Status ID	0032.000A	CS	U

Confidentiality Code	0040.1008	LO	U
Procedure Code Sequence	0008,1032	SQ	U
Reason For Performed Procedure Code Sequence	0040.1012	SQ	U
Other patient IDs Sequence	0010.1002	SQ	U
Institutional Department Name	0008.1040	LO	U
Number of Study Related Series	0020,1206	IS	U
Number of Study Related Instances	0020,1208	IS	U
Instance Availability	0008.0056	CS	U
SOP Classes in Study	0008,0062	UI	U
Current Patient Location	0038.0300	LO	U
Source Application Entity Title	0002.0016	AE	U
Time Zone Offset From UTC	0008.0201	SH	S

Types of Matching:

The types of Matching supported by the C-FIND SCP. An "S" indicates the identifier attribute uses Single Value Matching, an "R" indicates Range Matching, a "W" indicates wildcard matching, a 'U' indicates Universal Matching (i.e. return keys), and an 'L' indicates that UID lists are supported for matching. "NONE" means that single, range and wildcard matches are not supported on the key in question.

2.2.5.2 QIDO-RS Search for Series

Parameter	Restrictions
Media Types	The acceptable media types for the response payload. Examples include: "application/dicom+json" (send information in JSON format) "multipart/related; application/dicom+xml" (send information in XML format)
Matching Attributes	Refer to QIDO-RS Series Attribute Matching.
Relational Queries Supported	Yes
Hierarchical Search	Yes
Limit and Offset Supported	Yes

QIDO-RS Series Attribute Matching

Series Level - Attribute Name	Tag	VR	Types of Matching
Modality	0008,0060	CS	U
Series Description	0008,103E	LO	U
Series Instance UID	0020,000E	UI	U,S,W
Series Number	0020,0011	IS	U,S
Request Attributes Sequence > Requested Procedure ID	0040.0275 0040.1001	SQ SH	U,S
Request Attributes Sequence > Scheduled Procedure Step ID	0040.0275 0040.0009	SQ SH	U,S

Types of Matching:

The types of Matching supported by the C-FIND SCP. An "S" indicates the identifier attribute uses Single Value Matching, an "R" indicates Range Matching, a "W" indicates wildcard matching, a 'U' indicates Universal Matching (i.e. return keys), and an 'L' indicates that UID lists are supported for matching. "NONE" means that single, range and wildcard matches are not supported on the key in question.

2.2.5.3 QIDO-RS Search for Instances

Parameter	Restrictions
Media Types	The acceptable media types for the response payload. Examples include: application/dicom+json (send information in JSON format) multipart/related; type="application/dicom+xml" (send information in XML format)
Matching Attributes	Refer to QIDO-RS Series Attribute Matching.
Hierarchical Search	Yes
Relational Queries Supported	Yes
Limit and Offset Supported	Yes

QIDO-RS Instance Attribute Matching

Instance Level - Attribute Name	Tag	VR	Types of Matching
SOP Class UID	0008,0016	UI	U,W, L

SOP Instance UID	0008.0018	LO	U,W, L
Instance Number	0020.0013	IS	U,S

Types of Matching:

The types of Matching supported by the C-FIND SCP. An "S" indicates the identifier attribute uses Single Value Matching, an "R" indicates Range Matching, a "W" indicates wildcard matching, a 'U' indicates Universal Matching (i.e. return keys), and an 'L' indicates that UID lists are supported for matching. "NONE" means that single, range and wildcard matches are not supported on the key in question.

2.2.5.4 QIDO-RS Connection Policies

2.2.5.4.1 General

All standards RS connection policies apply. There are no extensions for RS options.

2.2.5.4.2 Number of Connections

By default, no limit on max simultaneous HTTPS request and max limit can be configured.

2.2.5.4.3 Response Status

Vitrea Connection – QIDO Service shall provide a response message with the following status code indicating different responses.

Code	Description
200	The request was successful
204	The request was successful but there are no matches
400	Bad request: the request is malformed. The query returns more results than the query limit set in the Image Archive service (if a query limit other than 0 (unlimited) has been set).
406	The Target Resource does not have a representation that would be acceptable to the user agent.
415	Invalid media type.
500	Internal service error.
503	Service unavailable.

2.2.5.4.4 Web Service Endpoint URL

<https://host:2525/rest/dicomweb>

2.2.6 Vitrea Connection WADO-URI Provider Specification

Vitrea Connection includes a WADO Server Application Entity that can be used to provide access to DICOM objects via a web based interface. This service may be configured to be solely accessible over TLS if required. The following table describes the configurable attributes of the WADO-URI Provider

WADO Configuration Parameters

Parameter	Configurable (Yes/No)	Default Value
Listening URL	Yes	http://localhost:8080/WADO

The following table describes the supported MIME types of the WADO-URI Provider.

WADO MIME TYPES

MIME Type	Description
application/dicom	The body content will be a “Part 10 File” that includes the meta-header. The returned DICOM object will be encoded using one the transfer syntaxes specified in the transfer syntax query parameter as defined below. By default, the transfer syntax shall be “Explicit VR Little Endian”.
image/jpeg	The body content will be a JPEG image. The size and quality of the image will depend on the query parameters as defined below.

The following table describes the query parameters that can be passed to the WADO-URI Provider.

WADO Query Parameters

Query Parameters	Required (Yes/No)	Description
requestType	Yes	The value must be “WADO”.
studyUID	Yes	The Study Instance UID identifying the DICOM object. The value must be present and encoded as a UID, except that it shall not be padded to even length with a NULL character.
seriesUID	Yes	The Series Instance UID identifying the DICOM object. The value must be present and encoded as a UID, except that it shall not be padded to even length with a NULL character.
objectUID	No	The SOP Instance UID identifying the DICOM object. If present, the value must be encoded as a UID, except that it shall not be padded to even length with a NULL character. If not present, the service will return a representative object (i.e. middle) from the series and study that are specified.
contentType	No	Must be image/jpeg.
rows	No	The maximum number of rows returned. Only valid for contentType image/jpeg.
columns	No	The maximum number of columns returned. Only valid for contentType image/jpeg.

2.2.7 Vitrea Connection WADO-RS Provider Specification

2.2.7.1 WADO-RS Retrieve Studies

Options	Restrictions
Data Types Supported (Accept Type)	multipart/related;type="application/dicom"
	multipart/related;type="application/octet-stream"
	multipart/related; type="image/dicom-rle" or the legacy multipart/related;type="image/dicom+rle"
	multipart/related;type = "image/jpeg" or the legacy multipart/related;type = "image/dicom+jpeg"
	multipart/related;type = "image/jls" or the legacy multipart/related;type = "image/dicom+jpeg-ls"
	multipart/related;type = "image/jp2" or the legacy multipart/related;type = "image/dicom+jp2"
	multipart/related;type = "image/jpx" or the legacy multipart/related;type = "image/dicom+jpx"
Transfer Syntaxes Supported (transfer-syntax Accept parameter)	Refer to the following: Transfer Syntaxes are represented by the term [Extended TS List]
SOP Class Restrictions	See SOP Classes for Vitrea Connection Storage Server AE

2.2.7.2 WADO-RS Retrieve Series

Options	Restrictions
Data Types Supported (Accept Type)	multipart/related;type="application/dicom"
	multipart/related;type="application/octet-stream"
	multipart/related; type="image/dicom-rle" or the legacy multipart/related;type="image/dicom+rle"
	multipart/related;type = "image/jpeg" or the legacy multipart/related;type = "image/dicom+jpeg"
	multipart/related;type = "image/jls" or the legacy multipart/related;type = "image/dicom+jpeg-ls"
	multipart/related;type = "image/jp2" or the legacy multipart/related;type = "image/dicom+jp2"
	multipart/related;type = "image/jpx" or the legacy multipart/related;type = "image/dicom+jpx"
Transfer Syntaxes Supported (transfer-syntax Accept parameter)	Refer to the following: Transfer Syntaxes are represented by the term [Extended TS List]
SOP Class Restrictions	See SOP Classes for Vitrea Connection Storage Server AE

2.2.7.3 WADO-RS Retrieve Instance

Options	Restrictions
Data Types Supported (Accept Type)	application/dicom
	multipart/related;type="application/dicom"
	multipart/related;type="application/octet-stream"
	image/dicom-rle or the legacy image/dicom+rle
	image/jpeg or the legacy image/dicom+jpeg
	image/jls or the legacy image/dicom+jpeg-ls

Options	Restrictions
	image/jp2 or the legacy image/dicom+jp2
	image/jpx or the legacy image/dicom+jpx
Transfer Syntaxes Supported (transfer-syntax Accept parameter)	Refer to the following: Transfer Syntaxes are represented by the term [Extended TS List]
SOP Class Restrictions	See SOP Classes for Vitrea Connection Storage Server AE

2.2.7.4 WADO-RS Retrieve Frames

Options	Restrictions
Data Types Supported (Accept Type)	multipart/related;type="application/octet-stream"
	multipart/related;type="image/dicom-rle" or the legacy multipart/related;type="image/dicom+rle"
	multipart/related;type="image/jpeg" or the legacy multipart/related;type="image/dicom+jpeg"
	multipart/related;type="image/jls" or the legacy multipart/related;type="image/dicom+jpeg-ls"
	multipart/related;type="image/jp2" or the legacy multipart/related;type="image/dicom+jp2"
	multipart/related;type="image/jpx" or the legacy multipart/related;type="image/dicom+jpx"
Transfer Syntaxes Supported (transfer-syntax Accept parameter)	Refer to the following: Transfer Syntaxes Supported (transfer-syntax Accept parameter)
SOP Class Restrictions	See SOP Classes for Vitrea Connection Storage Server AE

2.2.7.5 WADO-RS Retrieve Bulk Data

Options	Restrictions
Data Types Supported (Accept Type)	multipart/related;type="application/octet-stream"
	multipart/related;type="image/dicom-rle" or the legacy multipart/related;type="image/dicom+rle"
	multipart/related;type="image/jpeg" or the legacy multipart/related;type="image/dicom+jpeg"
	multipart/related;type="image/jls" or the legacy multipart/related;type="image/dicom+jpeg-ls"
	multipart/related;type="image/jp2" or the legacy multipart/related;type="image/dicom+jp2"
	multipart/related;type="image/jpx" or the legacy multipart/related;type="image/dicom+jpx"
Transfer Syntaxes Supported (transfer-syntax Accept parameter)	Refer to the following: Transfer Syntaxes Supported (transfer-syntax Accept parameter)
SOP Class Restrictions	See SOP Classes for Vitrea Connection Storage Server AE

2.2.7.6 WADO-RS Retrieve Metadata

Options	Restrictions
Data Types Supported (Accept Type)	multipart/related; type="application/dicom+xml"
	application/dicom+json
	application/vnd.canonmedical.dicom+fb
SOP Class Restrictions	See SOP Classes for Vitrea Connection Storage Server AE

2.2.7.7 WADO-RS Connection Policies

2.2.7.7.1 General

All standards RS connection policies apply. There are no extensions for RS options.

2.2.7.7.2 Number of Connections

By default, no limit on max simultaneous HTTPS request and max limit can be configured.

2.2.7.7.3 Web Service Endpoint URL

<https://host:2525/rest/dicomweb>

2.2.8 Vitrea Connection Proxy Server AE Specification

2.2.8.1 SOP Classes

Vitrea Connection Proxy Server Application Entity provides Standard Conformance to the following SOP Classes:

SOP Classes for Vitrea Connection Proxy Server AE

SOP Class Name	SOP Class UID	SCU	SCP
Verification SOP Class	1.2.840.10008.1.1	Yes	Yes
Patient Root Q/R Information Model – FIND	1.2.840.10008.5.1.4.1.2.1.1	Yes	Yes
Patient Root Q/R Information Model – MOVE	1.2.840.10008.5.1.4.1.2.1.2	Yes	Yes
Study Root Q/R Information Model – FIND	1.2.840.10008.5.1.4.1.2.2.1	Yes	Yes
Study Root Q/R Information Model – MOVE	1.2.840.10008.5.1.4.1.2.2.2	Yes	Yes
Modality Worklist Information Model – FIND	1.2.840.10008.5.1.4.31	Yes	Yes
Storage Commitment Push Model	1.2.840.10008.5.1.20.1	Yes	Yes
[Default SOP Class List]		Yes	Yes

2.2.8.2 Association Establishment Policy

2.2.8.2.1 General

The Vitrea Connection Proxy Server AE can both accept and propose Association Requests. The Vitrea Connection Storage Server AE will accept and propose Association Requests for the Verification, Storage, Query/Retrieve and Modality Worklist Query Services.

The DICOM standard application context name for DICOM 3.0 is always accepted and proposed:

DICOM application context name for Vitrea Connection Proxy Server AE

Application Context Name	1.2.840.10008.3.1.1.1
--------------------------	-----------------------

2.2.8.2.2 Number of Associations

Vitrea Connection Proxy Server can support multiple simultaneous Associations requested by peer Aes. Default is 10. This value can be configured through the configuration GUI.

Number of Associations accepted by Vitrea Connection Proxy Server AE

Maximum number of simultaneous Associations	10 (Configurable)
---	-------------------

2.2.8.2.3 Asynchronous Nature

Vitrea Connection Proxy Server does not support asynchronous communication. Multiple outstanding transactions are not supported. It allows up to one invoked and one performed operation on an Association (it is synchronous).

Asynchronous Nature as SCP for Vitrea Connection Proxy Server AE

Maximum number of outstanding asynchronous transactions	1 (Not Configurable)
---	----------------------

2.2.8.2.4 Implementation Identifying Information

The implementation information for this Application Entity is:

DICOM Implementation Class and Version for Vitrea Connection Proxy Server AE

Implementation Class UID	2.16.124.113638.2
Implementation Version Name	rialto-dcm4che2

2.2.8.3 Association Acceptance Policy

2.2.8.3.1 Activity – Receive Images

Description and Sequencing of Activities

A remote peer DICOM Application Entity, acting as an SCU, establishes an association with the Proxy Server AE that accepts these Associations for the purpose of receiving supported SOP Class Instances requests.

In the default configuration any Calling and Called AET will be accepted. If the Called AET does not correspond to the actual Proxy Server AET, only a Presentation Context for the Verification SOP Class will be accepted and the SCU can only verify the DICOM association but cannot invoke any other related DICOM service.

The Proxy Server AE may reject Association attempts as shown in the table below. The Result, Source and Reason / Diag columns represent the values returned in the corresponding fields of an ASSOCIATE-RJ PDU (see PS 3.8).

Association Rejection Reasons

Result	Source	Reason / Diag	Description
2 – rejected transient	provider	2 – local-limit exceeded	The (configurable) maximum number of simultaneous associations has been reached. An association request with the same parameters may succeed at a later time.
1 – rejected permanent	user	2 – application context name not supported	The Association request contained an unsupported Application Context Name. An association request with the same parameters will NOT succeed at a later time.
1 – rejected permanent	user	3 – calling AE title not recognized	The Association request contained an unrecognized Calling AE Title. An association request with the same parameters will NOT succeed at a later time unless configuration

changes are made.

Accepted Presentation Contexts

Vitreia Connection Proxy Server AE will accept Presentation Contexts as shown in the following table:

Accepted Presentation Contexts for Vitrea Connection Proxy Server AE

Abstract Syntax		Transfer Syntax	Role	Ext. Neg.
Name	UID			
Verification SOP Class	1.2.840.10008.1.1	[Default TS List]	SCP	None
Patient Root Q/R Information Model – FIND	1.2.840.10008.5.1.4.1.2.1.1	[Default TS List]	SCP	None
Study Root Q/R Information Model – FIND	1.2.840.10008.5.1.4.1.2.2.1	[Default TS List]	SCP	None
Patient Root Q/R Information Model – MOVE	1.2.840.10008.5.1.4.1.2.1.2	[Default TS List]	SCP	None
Study Root Q/R Information Model – MOVE	1.2.840.10008.5.1.4.1.2.2.2	[Default TS List]	SCP	None
Modality Worklist Information Model – FIND	1.2.840.10008.5.1.4.31	[Default TS List]	SCP	None
Storage Commitment Push Model	1.2.840.10008.1.20.1	[Default TS List]	SCP	None
[Default SOP Class List]			SCP	None

The preferred (default) Transfer Syntax is Explicit VR Little Endian.

SOP Specific Conformance

Vitreia Connection Proxy Server AE will forward the request to a single remote DICOM AE and return the response received from the remote DICOM AE. Therefore, all SOP specific conformances are inherited from the remote DICOM AE.

2.2.8.4 Association Initiation Policy

2.2.8.4.1 Activity – Synchronize Objects with an External Peer AE

Description and Sequencing of Activities

The Vitrea Connection Proxy Server AE will initiate a new Association when it receives an Association request from a remote DICOM AE.

Proposed Presentation Contexts

A list of Transfer Syntaxes can be configured for each AET. Vitrea Connection Proxy Server AE will propose Presentation Contexts as shown in the following table:

Abstract Syntax		Transfer Syntax	Role	Ext. Neg.
Name	UID			
Verification SOP Class	1.2.840.10008.1.1	[Default TS List]	SCU	None
Patient Root Q/R Information Model – FIND	1.2.840.10008.5.1.4.1.2.1.1	[Default TS List]	SCU	None
Patient Root Q/R Information Model – MOVE	1.2.840.10008.5.1.4.1.2.1.2	[Default TS List]	SCU	None

Study Root Q/R Information Model – FIND	1.2.840.10008.5.1.4.1.2.2.1	[Default TS List]	SCU	None
Study Root Q/R Information Model – FIND	1.2.840.10008.5.1.4.1.2.2.2	[Default TS List]	SCU	None
Modality Worklist Information Model – FIND	1.2.840.10008.5.1.4.31	[Default TS List]	SCU	None
Storage Commitment Push Model	1.2.840.10008.1.20.1	[Default TS List]	SCU	None
[Default SOP Class List]			SCU	None

SOP Specific Conformance

Vitreia Connection Proxy Server AE will forward the request to a single remote DICOM AE and return the response received from the remote DICOM AE. Therefore, all SOP specific conformance are inherited from the remote DICOM AE.

2.3 Network Interfaces

2.3.1 Supported Communication Stacks

DICOM Upper Layer over TCP/IP is supported.

2.3.1.1 TCP/IP Stack

Vitreia Connection inherits its TCP/IP stack from the installed Java Runtime Environment.

2.3.2 Physical Network Interface

Vitreia Connection is indifferent to the physical medium over which TCP/IP executes; it inherits this from the Java Runtime Environment.

2.4 Configuration

2.4.1 AE Title/Presentation Address Mapping

2.4.1.1 Local AE Titles

The local AE Titles and TCP ports are configurable through the Vitrea Connection configuration interface.

AE Title Configuration Table

Application Entity	Default AE Title	Default TCP/IP Port
Vitreia Connection Storage Server	RIALTO	11114
Vitreia Connection Query/Retrieve Server	RIALTO	11112
Vitreia Connection Foreign Exam Management Server	RIALTO	11111
Vitreia Connection STOW-RS Server	NA	8080
Vitreia Connection WADO URI Provider	NA	8080
Vitreia Connection Connect Server	RIALTO	4109
Vitreia Connection MWL Server	RIALTO	4106
Vitreia Connection Proxy Server	RIALTO	11116

Vitrea Connection DICOMweb (QIDO-RS and WADO-RS)	NA	2525
--	----	------

2.4.1.2 Remote AE Titles

Remote AE Titles, TCP/IP Addresses and ports can be configured through the Vitrea Connection configuration interface. In the default configuration, Association Requests with any Calling AET will be accepted.

2.4.2 Configuration Parameters

The following table shows the Vitrea Connection configuration parameters relevant to DICOM communication. Refer to the Vitrea Connection Service Manual for details on general configuration capabilities.

Configuration Parameter Table

Parameter	Configurable (Yes/No)	Default Value
General		
Maximum number of simultaneous Associations	Yes	10
Time-out waiting for A-ASSOCIATE RQ on open TCP/IP connection (ARTIM timeout)	Yes	5s
Time-out waiting on an open association for the next message (DIMSE timeout)	Yes	1h
Time-out waiting for acceptance or rejection Response to an Association Open Request. (Application Level timeout)	No	no timeout
Time-out waiting on an open association for the next message after sending A-RELEASE RSP or A-ABORT RQ (Closing timeout)	Yes	500ms
Maximum PDU size the AE can receive	Yes	16352
Maximum PDU size the AE can send	No	65535
Pack Command and Data PDVs in one PDU	Yes	false
Support for the Basic TLS Secure Transport Connection Profile	Yes	off
Accepted TLS Ciphers	Yes	-

3 Media Interchange

Vitreia Connection does not support Media Storage.

4 Support of Extended Character Sets

Vitreia Connection supports ISO_IR 100 (ISO 8859-1 Latin 1) as an extended character set.

5 Security

5.1 Security Profiles

Vitreia Connection supports secure DICOM communication in conformance with the Basic TLS Secure Transport Connection Profile. At default configuration, the TLS option is deactivated. In keeping with current security standards, TLS 1.0 and 1.1 are not supported.

When TLS is activated, the following ciphers are enabled by default:

- TLS_AES_128_GCM_SHA256
- TLS_AES_256_GCM_SHA384
- TLS_CHACHA20_POLY1305_SHA256
- ECDHE-ECDSA-AES128-GCM-SHA256
- ECDHE-RSA-AES128-GCM-SHA256
- ECDHE-ECDSA-AES256-GCM-SHA384
- ECDHE-RSA-AES256-GCM-SHA384
- ECDHE-ECDSA-CHACHA20-POLY1305
- ECDHE-RSA-CHACHA20-POLY1305
- DHE-RSA-AES128-GCM-SHA256
- DHE-RSA-AES256-GCM-SHA384

5.2 Association Level Security

Vitreia Connection can be configured to accept Association Requests from only a limited list of Calling AE Titles.

In the default configuration, Association requests with any Calling AET and any Called AET will be accepted. However, if the Called AET is not corresponding to any of the actual DICOM server AETs, only acceptance of the Presentation Context for Verification SOP Class will be returned in the Association Acceptance Response (A-ASSOCIATE AC).

5.3 Application Level Security

The Vitreia Connection web based administration module can be configured to require user authentication in order to access the user interface functionalities.